
mimoCoRB

Release 1.0.1

C. Mayer, K. Heitlinger, G. Quast

May 30, 2024

CONTENTS:

- 1 mimoCoRB - multiple-in multile-out Configurable Ring Buffer: Overview 1**
 - 1.1 For Developers: Description of components 1
 - 1.2 Simple application example 2
- 2 For application developers: Access Classes in the module *buffer_control* 7**
 - 2.1 Application examples 10
- 3 Module Documentation 19**
- Python Module Index 31**
- Index 33**

MIMOCORB - MULTIPLE-IN MULTILE-OUT CONFIGURABLE RING BUFFER: OVERVIEW

mimoCoRB: multiple-in multiple-out Configurable Ring Buffer

The package **mimoCoRB** provides a central component of each data acquisition system needed to record and pre-analyse data from randomly occurring processes. Typical examples are waveform data as provided by detectors common in quantum mechanical measurements, or in nuclear physics, particle physics and astro particle physics, e. g. photo tubes, Geiger counters, avalanche photo-diodes or modern SiPMs.

The random nature of such processes and the need to keep read-out dead times low requires an input buffer for fast collection of data and an efficient buffer manager delivering a constant data stream to the subsequent processing steps. While a data source feeds data into the buffer, consumer processes receive the data to filter, reduce, analyze or simply visualize the recorded data. In order to optimally use the available resources, multi-core and multi-processing techniques must be applied. Data consumers may be obligatory ones, i. e. data acquisition pauses if all input buffers are full and an obligatory consumer is still busy processing. A second type of consumers (random consumers or “observers”) receive an event copy from the buffer manager upon request, without interrupting the data acquisition process. Typical examples of random consumers are displays of a subset of the wave forms or of intermediate analysis results.

This project originated from an effort to structure and generalize data acquisition for several experiments in advanced physics laboratory courses at Karlsruhe Institute of Technology (KIT) and has been extensively tested under Ubuntu Linux.

As a simple demonstration, we provide data from simulated signals as would be recorded by a detector for cosmic muons with four detection layers. Occasionally, such muons stop in an absorber between the 2nd and 3rd layer, where they decay at rest and emit a high-energetic electron recorded as a 2nd pulse in one or two of the detection layers. After data acquisition, a search for typical pulse shapes is performed and data with detected double pulses are selected and copied into a second buffer. A third buffer receives data in a reduced format which only contains the parameters of accepted pulses. These data and the wave forms of all double-pulses are finally stored on disk. Such an application is a very typical example of the general process of on-line data processing in modern physics experiments and may serve as a starting point for own projects.

1.1 For Developers: Description of components

The following paragraphs provide some insight into the inner working of the components of *mimoCoRB* for interested users and for developers wanting to help improving the package. Application developers should use the more convenient access classes to build an application based the *mimoCoRB* framework.

In order to decouple the random occurrence of “events” one needs a buffer capable of rapidly storing new incoming data and delivering a constant data stream to subsequent consumer processes. This is typically implemented as a first-in, first out ringbuffer providing storage space in memory for incoming data, which is released and overwritten by new data when all consuming processes have finished.

As digital filtering of incoming data may be very CPU intensive, multi-processing and multi-core capable components are needed to ensure sufficient compute power to process and analyze data. *mimoCoRB.mimo_buffer* implements such a buffer allowing multiple processes to read (“multiple out”) or write (“multiple in”) to a shared buffer space.

Because processing of the data, i.e. digital filtering, selection, compression and storage or real-time visualization of the data can be a complex workflow, buffers may be arranged in chains where one or several reader processes of a buffer write to one or several output buffer(s).

The central component takes care of memory management and access control provided by the class **newBuffer**. To control the data flow in a full data acquisition suite, three types of access are needed, implemented as **Writer**, **Reader** and **Observer** classes. Readers of the same type are grouped together for multi-processing of compute-intense tasks and form a Reader-group. Observers receive only a sub-set of the data and are mainly intended to be used for visual inspection or graphical representation of samples of the recorded or processed data.

Processes for data provisioning from front-end hardware or from other sources, like disk files, web streams or simulation, rely on the Writer class; similarly, processes reading data for filtering and transfer to subsequent buffers or to analyse, transform or extract data to storage media use the Reader class. Note that the buffer manager ensures that every slot assigned to a Reader (or a group of Readers) is actually processed; therefore, input to a buffer blocks if the buffer is filled up completely. The Writer class resumes data input as soon as a Reader or member of a Reader-group has finished processing and thus freed a slot in the buffer.

Multiprocessing is enabled by use of the *shared_memory* module of the *multiprocessing* package available since Python 3.8 for direct access to shared memory across processes. Other modules of the package (*Process*, *Lock*, *Event*, and *SimpleQueue* or *Queue*) are used to create and control sub-processes and for signalling and message or data exchange across processes.

The format of data stored in the buffers is based on structured *numpy* arrays with (configurable) field names and *numpy* *dtypes*. Each buffer entry is also associated with a unique number, a time stamp and a deadtime fraction to be provided by the initial data producer. The deadtime accounts for inefficiencies of the data acquisition due to processing in *mimoCoRB*. These metadata are set by the initial producer and must not be changed at a later stage in the processing chain.

1.2 Simple application example

An application example of *mimo_buffer* is shown below; it is also provided as a unit test. The set-up is as follows:

Two ring buffers are defined:

- input Buffer RB_1: 10 ch x 1024 slots (int32)
- output Buffer RB_2: 10 ch x 2 slots/ch (float64)

Simple data is filled into RB_1, copied and extended by a process writing data into RB_2, and finally a reader process to check integrity and completeness of the data. The most complex part of the code is in function *run_control()*, which demonstrates how to set up the buffers, define Reader and Writer instances and start the parallel processes for generating, processing and reading the data.

The example including comment lines for explanation is shown here:

```
import time
import unittest
import numpy as np
from multiprocessing import Process, Value
from mimocorb import mimo_buffer as bm

# global variables
N_requested = 1000 # number of data injections ("events")
```

(continues on next page)

(continued from previous page)

```

Time_tick = 0.001    # time between events
Ncpul = 2            # number of parallel analyzer processes

def data_generator(sink_dict):
    """writes continuously rising integers to buffer specified in sink_dict
    """
    sink = bm.Writer(sink_dict)
    n=0
    # inject data
    for x in range(N_requested):
        buffer = sink.get_new_buffer() # get new buffer and pass last item
        # random wait for next data item
        time.sleep(-Time_tick*np.log(np.random.rand() ))
        # fill "data"
        n += 1
        buffer[:] = n
    # process last data item
    sink.process_buffer()

def analyzer(source_dict, sink_dict):
    """read from source and write first element and a time difference to sink
    """
    source = bm.Reader(source_dict)
    sink = bm.Writer(sink_dict)
    start_time = time.time()

    while True:
        input_data = source.get()
        output_data = sink.get_new_buffer()
        # process data
        output_data[0] = input_data[0]
        # mimick processing time
        time.sleep(2*Time_tick)
        output_data[1] = time.time() - start_time

        #
        sink.process_buffer()

def check_result(source_dict, res):
    """reads RB_2 and sum up the integer content

    sum is returned as shared memory Value-object
    """
    source = bm.Reader(source_dict)
    sum_rb = 0
    while True:
        input_data = source.get()
        res.value +=int(input_data[0])

def run_control():
    """Setup buffers, start processes and shut_down when 1st writer done
    """

```

(continues on next page)

(continued from previous page)

```

# Create ring buffers: #2: 10 channel, 2 value per channel
#   (1: buffer content; 2: time difference as int)
#   d_type = [('chA', np.float)] #not necessary: always the same type
generator_buffer = bm.NewBuffer(10, 1, np.int32)
eval_buffer = bm.NewBuffer(10, 2, np.float32)

# create readers first
source_dic_gen = generator_buffer.new_reader_group()
source_dic_eval = eval_buffer.new_reader_group()

# Create worker processes (correct sequence: first action as last)
process_list = []
# evaluation to test ring buffer behavior
result = Value('i', 0) # int variable in shared memory
process_list.append(Process(target=check_result,
                           args=(source_dic_eval, result)))
# data transfer between the 2 buffers: generator_buffer -> eval_buffer
sink_dic_eval = eval_buffer.new_writer()
# work with all cpu's requested
number_of_workers = Ncpu1
for i in range(number_of_workers):
    process_list.append(Process(target=analyzer,
                              args=(source_dic_gen, sink_dic_eval)))

# fill buffer (generator_buffer) with data first
sink_dic_gen = generator_buffer.new_writer()
process_list.append(Process(target=data_generator,
                           args=(sink_dic_gen,)))

for p in process_list:
    p.start()

run_active = True
while run_active:
    run_active = False if process_list[-1].exitcode==0 else True
    time.sleep(0.1) # wait
time.sleep(0.1) # some grace-time for readers to finish

generator_buffer.shutdown()
eval_buffer.shutdown()
del generator_buffer, eval_buffer

for p in process_list:
    p.join()

return result.value

class RPTest(unittest.TestCase):

    def test_process(self):
        # start python test module and check result

```

(continues on next page)

(continued from previous page)

```
a = run_control()
expected_result = N_requested*(N_requested+1)//2
    # expected result: sum(i); i = 1, N_requested
self.assertEqual(a, expected_result)

if __name__ == "__main__":
    unittest.main(verbosity=2)
```

FOR APPLICATION DEVELOPERS: ACCESS CLASSES IN THE MODULE *BUFFER_CONTROL*

To facilitate user interaction with the buffer manager a set of additional classes is provided in the module *buffer_control* to set-up and manage cascades of ringbuffers and the associated functions for filling, filtering and extracting data. These classes are most interesting for application developers wanting to build upon the *mimoCoRB* framework.

The classes are:


- **class *buffer_control***
Set-up and management of ringbuffers and associated sub-processes. This is the overarching class with access to all created buffers and sub-processes.
- **class *rbImport***
Read data from source (front-end like a PicoScope USB oscilloscope, or from file or simulation) and put data in a *mimo_buffer*. Data input is handled by a call of a user-supplied Python generator (i.e. via 'yield()') for data and metadata.
- **class *rbTransfer***
Read data from a *mimo_buffer*, filter and/or reformat data and write to output *mimo_buffer(s)*. Data is provided as the argument to a user-defined filter function returning *None* if data is to be discarded, a number if data is to be copied to another buffer, or - optionally - a list of transformed data records produced from processed input data. If such data are provided, a respective number of ringbuffers as destination must be configured.
- **class *rbExport***
Read data from *mimo_buffer* and analyze (with user-supplied code), without writing to another ringbuffer. Data is expected to be provided by a Python generator in the `__call__()` method of the class yielding a tuple of data and metadata.
- **class *rbObserver***
Deliver data from a buffer to an observer process. A tuple (data, metadata) is provided by a Python generator implemented in the `__call__()` method of the class.
- **class *rb_toTxtfile*:**
Save *mimo_buffer* data to a file in csv-format. The header line of this file contains the keys of the respective columns, which are derived from the datatype of the structured ringbuffer array. Aliases for improved clarity can be provided in the configuration file.
- **class *rb_toParquetfile*:**
Save *mimo_buffer* data to an archive in *tar* format; each data record is packed in Parquet format.
- **class *run_mimoDAQ***
Setup and run a Data Acquisition suite with the *mimoCoRB* buffer manager. The layout of ringbuffers and associated functions are defined in a configuration file in *yaml* format. All configured functions are executed as worker processes in separate sub-processes and therefore optimal use is made of multi-core architectures.

- *class **bufferinfoGUI**:*

A graphical interface showing buffer rates and status information and providing some control buttons interacting with the `run_mimoDAQ` class.


These classes shield much of the complexity from the user, who can thus concentrate on writing the pieces of code need to acquire and process the data. The access classes expect as input lists of dictionaries with the parameters of buffers to read from (**source_list**), to write to (**sink_list**) or to observe (**observe_list**). An additional dictionary (**config_dict**) provides the parameters needed for the specific functionality, for example names of functions to read, filter or manipulate data or the names of target files. The interface for passing data between the user-defined functions and ringbuffers relies on Python generators (i.e. the `yield` instruction).

The overarching class **buffer_control** provides methods to setup buffers and worker processes and to control the data acquisition process. The methods collected in the class `run_mimoDAQ`, in particular the function `run_mimoDAQ`, contains the code needed to run a real example of a data-acquisition suite defined in a configuration file specifying the associated, user-defined functions for data provisioning, filtering and storage. `run_mimoDAQ` is controlled either by keyboard commands or from a graphical user interface; pre-defined conditions on the total number of events processed, the duration of the data taking run or finishing of the writer process to the first buffer due to source exhaustion can also be defined to end data taking. The class structure and dependencies are shown in the figure below.

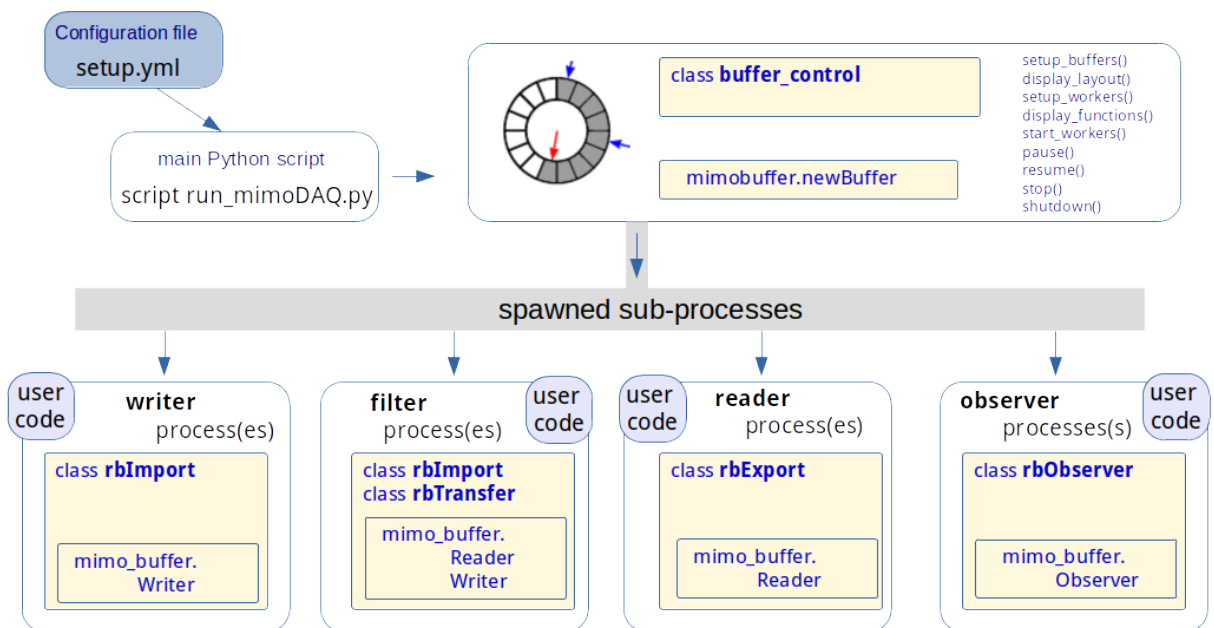

A sequence diagram of a typical application, shown below, illustrates the interplay and dependencies of the classes described above. The script `run_daq.py` creates an instance of `run_mimoDAQ` and starts its `run()`-method. The interaction with the user-supplied filter functions is handled by methods of the class `buffer_control`.

mimocorb sequence (main parts)

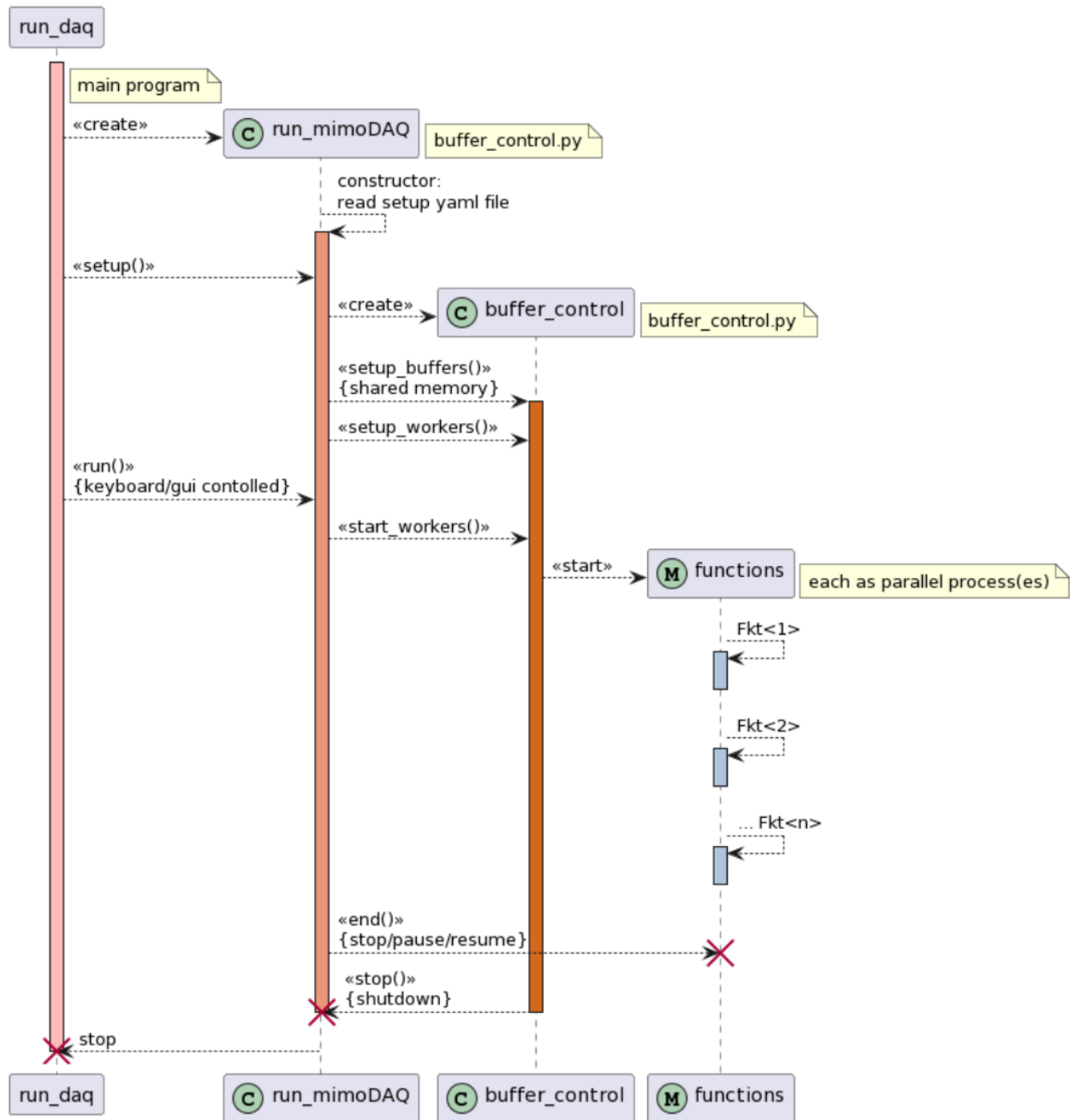

For complex setups and longer data-taking periods it is important to gain a quick overview of the status of all buffers and to monitor long-term stability. Therefore, a graphical display with the processing rate of all buffers is provided by the class **bufferinfoGUI**. A text window receives frequent updates of the number of events processed by each buffer and of the buffer fill-levels. Clickable control buttons send information via a dedicated command queue to the calling process *run_mimoDAQ* and enable pausing, resuming and controlled ending of the data-acquisition processes.

The suggested structure of the project work-space for mimoCoRB applications is as follows:

```
|--> <user working directory>          # the main configuration script resides here
```

(continues on next page)

(continued from previous page)

```
| --> modules      # project-specific, user-supplied python code
| --> config       # configuration files in yaml format
| --> target       # output of data-acquisition run(s)
```

For illustration and as a starting point for own applications, a stand-alone example is provided as part of the package, as described in the following section.

2.1 Application examples

The subdirectory *examples/* contains some complete application use cases based on input waveforms that are simulated in real-time.

The central piece of every *mimoCoRB* application is the configuration file; examples of different complexity are provided in *examples/*_setup.yaml*. Code snippets for data input, filtering and output and as well as configuration files are provided in the subdirectories *examples/modules/* and *examples/config/*, respectively.

Waveform data, as provided by, for example, a multi-channel digital oscilloscope, are generated and filled into the first one of a cascaded set of ringbuffers and passed on to subsequent buffer stages depending on filter conditions. Data in the last buffer are recorded to disk. The configuration files and the recorded data files are stored in the subdirectory *examples/target/<projectname>_<date_and_time>*.

All examples run stand-alone and use as input simulated waveform data of short pulses as they arise e.g. in scintillator detectors. The simulated physics process corresponds to signatures produced by cosmic muons penetrating the detector layers.

Simple Example

A very simple example consists of recording two input channels from two redundant detectors (e.g. two stacked muon panels). The buffer configuration is defined in the file *examples/spectrum_setup.yaml*, shown here:

```
# Configuration for recording two channels with mimoCoRB
# -----
#
# configure two buffers:
# - RB_1 for raw waveforms
# - RB_2 for derived pulse parameters
# data from RB_2, the result buffer, are saved to a file in csv (text) format.
#
# - data from RB_1 are also passed to an Observer process driving a real-time display
# - data RB_2 are passed to a Reader process driving a real-time histogram display
#
# Notes:
#
# 1. additional config files controlling the user functions are
#    located in the subdirectory config/
# 2. necessary application-specific user code is located
#    in the subdirectory modules/
# -----
#
RingBuffer:
# define ring buffers
```

(continues on next page)

(continued from previous page)

```

- RB_1:
  # raw input data buffer (from picoScope, file or simulation)
  number_of_slots: 16
  channel_per_slot: 500
  data_type:
    1: ['chA', "float32"]
    2: ['chB', "float32"]
- RB_2:
  # buffer with correct signature double pulse parameters
  number_of_slots: 16
  channel_per_slot: 1
  data_type:
  data_type:
    1: ['chA_height', "float32"]
    2: ['chA_position', "int32"]
    3: ['chA_integral', "float32"]
    4: ['chB_height', "float32"]
    5: ['chB_position', "int32"]
    6: ['chB_integral', "float32"]

```

Functions:

```

# define functions and ringbuffer assignment
- Fkt_main:
  config_file: "config/spectrum_config.yaml"
- Fkt_1:
  ## for simulation source
  file_name: "modules/simul_source"
  fkt_name: "simul_source"
  num_process: 1
  RB_assign:
    RB_1: "write"
- Fkt_2:
  file_name: "modules/spectrum_filter"
  fkt_name: "find_peaks"
  num_process: 2
  RB_assign:
    RB_1: "read"
    RB_2: "write"
- Fkt_3:
  file_name: "modules/save_files"
  fkt_name: "save_to_txt"
  num_process: 1
  RB_assign:
    RB_2: "read"

```

The example coming with this package contains two more convenience functions, one for an observer process displaying a random sample of waveforms in an oscilloscope display, and a second one for on-line analysis and histogramming of buffer data. The addendum to the configuration looks as follows:

```

# --- the following functions are optioal
- Fkt_4:
  file_name: "modules/plot_waveform"
  fkt_name: "plot_waveform"

```

(continues on next page)

(continued from previous page)

```

num_process: 1
RB_assign:
  RB_1: "observe"
- Fkt_5:
  file_name: "modules/plot_histograms"
  fkt_name: "plot_histograms"
  num_process: 1
  RB_assign:
    RB_2: "read" # pulse parameters

```

These additional functions rely on the modules *mimocorb.plot_buffer* and *mimocorb.histogram_buffer*, which provide animated displays of waveforms similar to an oscilloscope and a histogram package for life-updates of frequency distributions of scalar variables.

Configuration parameters needed for the functions associated to the ringbuffers can either be specified as a *yaml* block under a keyword in the general configuration file that is assigned to the function *Fkt_main* in the example above (*config_file*: "*config/spectrum_config.yaml*").

The functions are started as sub-processes and have a unique interface. Lists of dictionaries provide the necessary information to connect to the buffer manager via the *Writer*, *Reader* or *Observer* classes of the package. This information comprises the pointer to the shared buffer manager as well as pointers to instances of the functions *Event()* or *Queue()* from the multiprocessing package to enable communication and data transfers across processes. A further dictionary (*config_dict*) provides the function-specific parameters discussed previously. The keyword dictionary *rb_info* specifies whether writer, reader or observer functionality is required. It contains a copy of the ring-buffer assignment block („RB_assign:”) from the main setup file on function level. Its purpose is to facilitate the ring-buffer access part within a function or class.

The function interface looks as follows:

```

def <function_name>(  

    source_list=None, sink_list=None, observe_list=None, config_dict=None, **rb_info):

```

This interface must be respected by any user function. The argument list must also be passed to instances of the access classes *rbImort*, *rbExport*, *rbTransfer* or *rbObserve*. An example of a user function in the directory *modules/* to write buffer data to a text file is shown below:

```

"""Module save_files to handle file I/O for data in txt and parquet format

    This module relies on classes in mimocorb.buffer_control
"""
from mimocorb.buffer_control import rb_toTxtfile, rb_toParquetfile
def save_to_txt(source_list=None, sink_list=None, observe_list=None, config_dict=None,
    ↪ **rb_info):
    sv = rb_toTxtfile(source_list=source_list, config_dict=config_dict, **rb_info)
    sv()

```

Starting this example with the command *../run_daq.py spectrum_setup.yaml* yields the following output on screen:

```

*** script ../run_daq.py running

2 buffers created... List of buffers
RB_1 16 500
RB_2 16 1
List of functions

```

(continues on next page)

(continued from previous page)

```

FKT_1 simul_source (1)    {'RB_1': 'write'}
FKT_2 find_peaks (2)     {'RB_1': 'read', 'RB_2': 'write'}
FKT_3 save_to_txt (1)    {'RB_2': 'read'}
FKT_4 plot_waveform (1)  {'RB_1': 'observe'}
FKT_5 plot_histograms (1) {'RB_2': 'read'}

```

Two buffers are created in this case, *RB_1* and *RB_2*. *RB_1* with two channels with 500 samples each is the input buffer, *RB_2* contains 6 scalar variables and is the output buffer, which is filled by the function *find_peaks* with two active workers. The functions *save_to_txt* and *plot_histograms* read from this buffer and store data to dist or show histograms, respectively. The function *plot_waveform* takes random samples from *RB_1* and displays the raw waveform data.

This example serves as a convenient starting point for own application development. The code in *simul_source.py*, shown below, is a very general example for data input. Only the function *pulseSimulator()* needs to be replaced by a function providing data from your own source.

```

"""**simul_source**: a simple template for a mimoCoRB source to
enter simulated waveform data in a mimoCoRB buffer.

Input data is provided as numpy-array of shape (number_of_channels, number_of_samples).
"""

from mimocorb.buffer_control import rbImport
import numpy as np
import sys, time
from pulseSimulator import pulseSimulator

def simul_source(source_list=None, sink_list=None, observe_list=None, config_dict=None,
↳ **rb_info):
    """
    Generate simulated data and pass data to buffer
    The class mimocorb.buffer_control/rbImport is used to interface to the
    newBuffer and Writer classes of the package mimoCoRB.mimo_buffer

    This example may serve as a template for other data sources

    :param config_dict: configuration dictionary

    - events_required: number of events to be simulated or 0 for infinite
    - sleeptime: (mean) time between events
    - random: random time between events according to a Poission process
    - number_of_samples, sample_time_ns, pretrigger_samples and analogue_offset
    describe the waveform data to be generated (as for oscilloscope setup)
    """

    events_required = 1000 if "eventcount" not in config_dict else config_dict[
↳ "eventcount"]

    def yield_data():
        """generate simulated data, called by instance of class mimoCoRB.rbImport"""

        event_count = 0
        while events_required == 0 or event_count < events_required:
            pulse = dataSource(number_of_channels)

```

(continues on next page)

(continued from previous page)

```

# deliver pulse data and no metadata
yield (pulse, None)
event_count += 1

dataSource = pulseSimulator(config_dict)
simulsource = rbImport(config_dict=config_dict, sink_list=sink_list, ufunc=yield_
↳data, **rb_info)
number_of_channels = len(simulsource.sink.dtype)

simulsource()


```

Complex Example

In the second, more complex example discussed here we consider multiple pulses on the same channels, where, e.g. the first pulse originates from a detected muon and the second, later one from a decay electron of a muon that has been stopped in or near a detection layer.

The raw data are analyzed, and accepted data with a double-pulse signature are selected and directly passed on to a second ringbuffer. A third buffer contains only the information on found signal pulses; a result file in *csv* format contains the data extracted from this buffer.

A graphical representation of the set-up is shown in the figure below [source: Master's Thesis Christoph Mayer, ETP 2022]. Note that the oscilloscope is replaced by a signal simulation in the example provided.

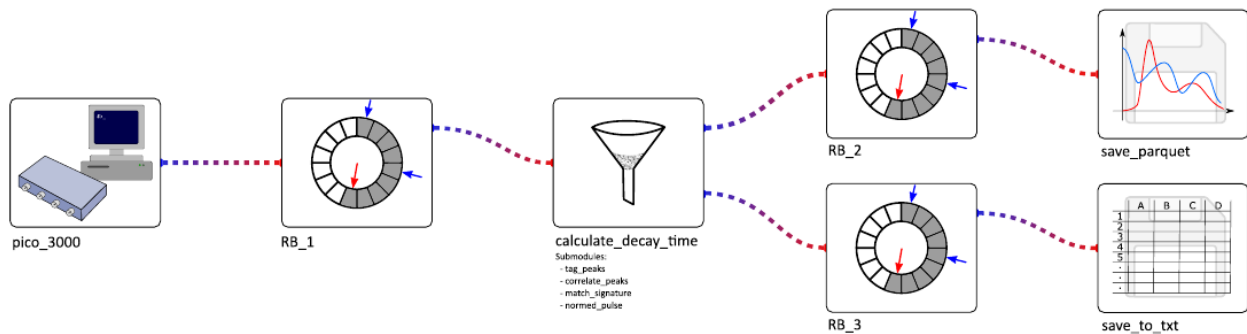

The buffer layout and the associated functions are defined in the main configuration file *simulsource_setup.py*, which serves as the input to the execution script *run_daq.py* in the top-level directory of the package. The *python* files *simul_source.py*, *lifetime_filter.py* and *save_files.py* contain the user code for data generation, analysis and filtering and extraction of the finally accepted data to disk files. The *.yaml* files *simulation_config.yaml* and *save_lifetimes.yaml* contain configurable parameters provided to these functions.

This example is executed from the directory *examples/* by entering:

```
../run_daq.py simulsource_setup.yaml
```

Again, the screen output gives an overview of the generated buffers and the functions writing to and reading from them:

```

*** script ../run_daq.py running

3 buffers created... List of buffers
RB_1 128 4250
RB_2 128 4250
RB_3 32 1
List of functions

```

(continues on next page)

(continued from previous page)

```

FKT_1 simul_source (1)  {'RB_1': 'write'}
FKT_2 calculate_decay_time (2)  {'RB_1': 'read', 'RB_2': 'write', 'RB_3': 'write'}
FKT_3 save_to_txt (1)  {'RB_3': 'read'}
FKT_4 save_parquet (1)  {'RB_2': 'read'}
FKT_5 plot_waveform (1)  {'RB_2': 'observe'}
FKT_6 plot_histograms (1)  {'RB_3': 'read'}

```

The input *yaml* file for this example looks as follows:

```

# Application example for mimoCoRB
# -----
#
# three buffers:
# - RB_1 for (simulated) raw waveforms
# - RB_2 for selected double-pulses
# - RB_3 for derived pulse parameters
#
# data from RB_2 and RB_3 are saved to files in tarred parquet format
# or in text format.
#
# data from RB_2 are passed to an observer process driving a real-time display
#
# Notes:
#
# 1. additional config files controlling the user functions are
#    located in the subdirectory config/
# 2. user necessary application-specific user code is located
#    in the subdirectory modules/
# -----
#
# general control options
output_directory: target      # directory to store output data
GUI_control: true             # control daq via Grapical User Interface
KBD_control: true             # control daq via KeyBoarD

RingBuffer:
  # define ring buffers
  - RB_1:
      # raw input data buffer (waveforms from PicoScope, file_source or simulation)
      number_of_slots: 128
      channel_per_slot: 4250
      data_type:
        1: ['chA', "float32"]
        2: ['chB', "float32"]
        3: ['chC', "float32"]
        4: ['chD', "float32"]
  - RB_2:
      # buffer with accepted signatures (here double-pulses)
      number_of_slots: 128
      channel_per_slot: 4250

```

(continues on next page)

(continued from previous page)

```

    data_type:
      1: ['chA', "float32"]
      2: ['chB', "float32"]
      3: ['chC', "float32"]
      4: ['chD', "float32"]
- RB_3:
  # buffer with pulse parameters (derived from waveforms)
  number_of_slots: 32
  channel_per_slot: 1
  data_type:
    1: ['decay_time', "int32"]
    3: ['1st_chA_h', "float32"]
    4: ['1st_chB_h', "float32"]
    5: ['1st_chC_h', "float32"]
    6: ['1st_chA_p', "int32"]
    7: ['1st_chB_p', "int32"]
    8: ['1st_chC_p', "int32"]
    9: ['1st_chA_int', "float32"]
    10: ['1st_chB_int', "float32"]
    11: ['1st_chC_int', "float32"]
    12: ['2nd_chA_h', "float32"]
    13: ['2nd_chB_h', "float32"]
    14: ['2nd_chC_h', "float32"]
    15: ['2nd_chA_p', "int32"]
    16: ['2nd_chB_p', "int32"]
    17: ['2nd_chC_p', "int32"]
    18: ['2nd_chA_int', "float32"]
    19: ['2nd_chB_int', "float32"]
    20: ['2nd_chC_int', "float32"]
    21: ['1st_chD_h', "float32"]
    22: ['1st_chD_p', "int32"]
    23: ['1st_chD_int', "float32"]
    24: ['2nd_chD_h', "float32"]
    25: ['2nd_chD_p', "int32"]
    26: ['2nd_chD_int', "float32"]
Functions:
  # define functions and assignments
- Fkt_main:
  config_file: "config/simulation_config.yaml"
- Fkt_1:
  file_name: "modules/simul_source"
  fkt_name: "simul_source"
  num_process: 1
  RB_assign:
    RB_1: "write"
- Fkt_2:
  file_name: "modules/lifetime_filter"
  fkt_name: "calculate_decay_time"
  num_process: 2
  RB_assign:
    RB_1: "read"      # input
    RB_2: "write"     # waveform to save (if double pulse was found)

```

(continues on next page)

(continued from previous page)

```

        RB_3: "write"      # pulse data
- Fkt_3:
  file_name: "modules/save_files"
  fkt_name: "save_to_txt"
  config_file: "config/save_lifetime.yaml"
  num_process: 1
  RB_assign:
    RB_3: "read"        # pulse data
- Fkt_4:
  file_name: "modules/save_files"
  fkt_name: "save_parquet"
  num_process: 1
  RB_assign:
    RB_2: "read"        # waveform to save

```

The functions for plotting result variables and for histogramming are also contained in the example configuration provided as part of the package, in the same way as explained for the first example.

In addition to the remarks concerning configuration files, it is worth noticing here that a dedicated configuration file is specified in a separate *yaml* file for function *Fkt_3*. This possibility is particularly useful if the same function code is used to handle data from different buffers, e.g. the file name and special formatting for writing buffer contents to a file in *csv* format.

Even more omplex Example

A similar, but even more complex case with two streams going to different output buffers and files is specified in the file *simul_spin_setup.yaml*. Starting this example with `> ./run_daq.py simul_spin_setup.yaml` leads to the following terminal output:

```


*** script ../run_daq.py running

4 buffers created... List of buffers
RB_1 128 4250
RB_2 128 4250
RB_3 32 1
RB_4 32 1
List of functions
FKT_1 simul_source (1) {'RB_1': 'write'}
FKT_2 calculate_decay_time (2) {'RB_1': 'read', 'RB_2': 'write', 'RB_3': 'write', 'RB_
↪4': 'write'}
FKT_3 save_to_txt (1) {'RB_3': 'read'}
FKT_4 save_to_txt (1) {'RB_4': 'read'}
FKT_5 save_parquet (1) {'RB_2': 'read'}
FKT_6 plot_waveform (1) {'RB_2': 'observe'}
FKT_7 plot_histograms (1) {'RB_3': 'read'}
FKT_8 plot_histograms (1) {'RB_4': 'read'}

Running 260s RB_1: 76409(0) 302Hz RB_2: 6118(0) 19Hz RB_3: 2952(1) 9.99Hz RB_4: 3166(1)
↪8.99Hz

```

A screenshot of a data-acquisition run with input from simulated data is shown in the figure below.

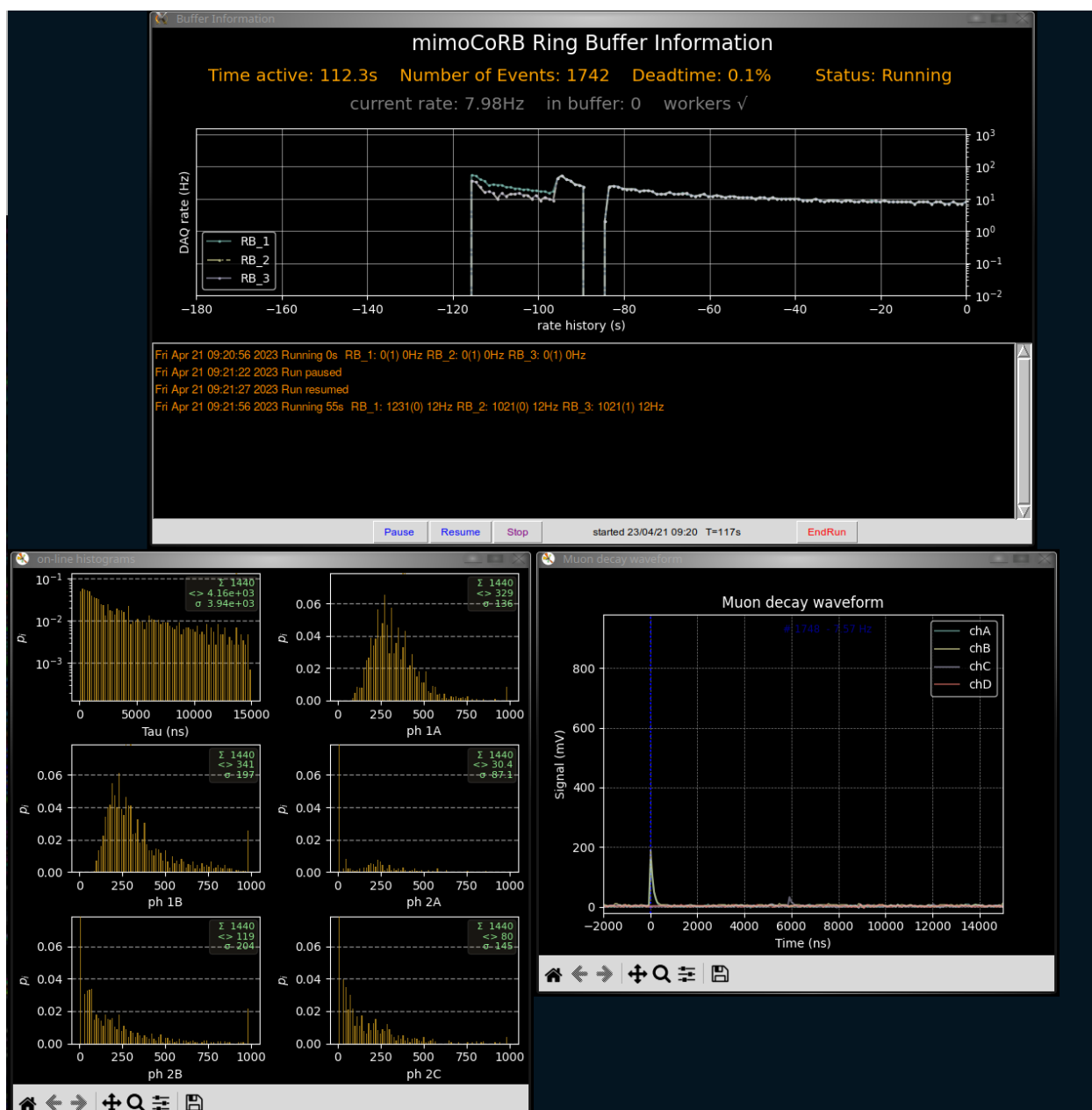

MODULE DOCUMENTATION

mimo-ringbuffer:

Module implementing a multiple-in multiple-out ringbuffer appropriate for multi-core multiprocessing.

The ringbuffer creation and management is handled by the class `NewBuffer`. Buffer access is managed by the `Reader`, `Writer` and `Observer` classes.

classes:

- `NewBuffer`: create a new ringbuffer, assign writer(s) and reader(s) or observer(s)

methods:

- `new_writer`
 - `new_reader_group`
 - `new_observer`
 - `buffer_staus`
- `Writer`: write elements into a ringbuffer
 - `Reader`: read all elements from a ringbuffer
 - `Observer`: read selected elements from a ringbuffer.

class `mimocorb.mimo_buffer.NewBuffer`(*number_of_slots, values_per_slot, dtype, debug=False*)

Class to create a new ringbuffer object according to the 'FIFO' principle (first-in first-out).

Memory shares, IPC queues, lock and event objects as well as background threads are defined for the multiprocessing ringbuffer management. Methods are provided to build the setup dictionaries (necessary parameter objects) for the `Reader`, `Writer` or `Observer` instances, respectively. Further, methods are provided to allow an index processing (e.g. listeners) and to pause data processing.

Index processing overview: to achieve a proper handling for writing and reading data into a ringbuffer slot regarding the possible definition of more than one reader for a ringbuffer (multiprocessing mode on function level) the determination of the correct slot index is essential. For a `Writer` or `Reader` two `SimpleQueues` are defined, respectively. In principle one queue contains the index to be used next and the other the index just processed. An `Observer` is treated in a simpler way by using an own global index variable for each instance.

The queues are passed via the setup dictionary to the corresponding `Reader` or `Writer` instance, respectively, and processed there. In addition, the buffer manager provides methods to control the index determination (`_writer/_reader/_observer_queue_listener()`), started in own threads.

Writer index:

- `writer_empty_queue`: contains the slot numbers (initially filled); defines the next free ringbuffer slot.
 - fetched (removed) in the class `Writer` -> `get_new_buffer()`

- last processed slot number is refilled in `_increment_reader_pointer()`
- `writer_filled_queue`: empty; contains the last processed slot number (distributed to all defined readers).
 - `process_buffer()`

Reader index:

- manually incremented in `_increment_reader_pointer()` via global variable `read_pointer`
- `done_queue`: empty; already processed slot number
- fetched in `_reader_queue_listener()`
- filled in the class `Reader` -> `get()` via the global variable `_last_get_index`
- `todo_queue`: empty; slot number to be processed next
- fetched in the class `Reader` -> `get()`
- filled in `_writer_queue_listener` (within the list `reader_todo_queue_list`)

Observer index:

- the global variable `obs_pointer` is used; it is an early copy of the `write_pointer` variable
 - defined in `_writer_queue_listener()`
 - directly used as index in `_observeQ_listener()`

important methods:

- `__init__()` constructor to create a new 'FIFO' ringbuffer
- `new_writer()` create new writer
- `new_reader_group()` create reader group
- `new_observer()` create observer
- `buffer_status()` display status: event count, processing rate, occupied slots
- `pause()` disable writer(s) to ringbuffer
- `resume()` (re-)enable writers
- `set_ending()` stop data-taking (gives processes time to finish before shutdown)
- `close()` release shared memory
- `shutdown()` end connected processes, delete ringbuffer

buffer_status()

Processing Rate and approximate number of free slots in this ringbuffer. This method is meant for user information purposes only, as the result may not be completely accurate due to race conditions.

Returns

cumulative event count, number of free slots, processing rate, average deadtime

Return type

tuple

new_observer()

Method to create a new (Queue based) observer.

Method: a copy of the most recent data (latest write_pointer) is transferred via a Queue whenever the Queue (of size 1) is empty. Sending data through the Queue is handled in a sparate thread

Returns

The `setup_dict` object passed to an `Observer`-instance to give access to the data Queue defined for this ringbuffer.

Return type

dict

new_reader_group()

Method to create a new reader group. The processing workload of a group can be distributed to multiple processes by using the same setup dictionary (`setup_dict`) defined for a `Reader`-object. Each ringbuffer element is processed by one reader group process. It's possible to create multiple reader groups per ringbuffer, where each reader group gets every element written to the ringbuffer. If a reader group is created, at least one `Reader`-class instance MUST steadily call its `get()` method to prevent the ringbuffer from blocking and to allow a safe shutdown.

Returns

The `setup_dict` object passed to a `Reader`-instance to grant read access to this ringbuffer.

Return type

dict

new_writer()

Method to create a new writer. It is possible to create multiple writers and simply share a setup dictionary definition between different `Writer`-instances (analogues to the behavior of the `new_reader_group`).

Returns

The `setup_dict` object passed to a `Writer`-instance to grant write access to this ringbuffer.

Return type

dict

pause()

Disable writing to ringbuffer (paused)

resume()

(Re)enable writing to ringbuffer (resume)

set_ending()

Stop data flow (before shut-down)

shutdown()

Shut down the ringbuffer(s): close background threads, terminate associated processes and release the shared memory definitions.

Affect processes using a `Reader`, `Writer` or `Observer` instance to a ringbuffer.

A 'trickle down' approach is used to have as few ringbuffer elements as possible unprocessed. This may not work correctly with more complex signal analysis chains. So always make sure to shut down the ringbuffers in data flow order (start with first element of the chain, the ringbuffer closest to the signal source).

CAUTION! If there are loops in the signal analysis chain, this method may end in an infinite loop!

class mimocorb.mimo_buffer.**Observer**(*setup_dict*)

Class for reading selected elements from a ringbuffer via q multiprocessing Queue

The data transfer is implemented via a multiprocessing Queue and interfaces with the ringbuffer manager (NewBuffer-class).

get()

Get latest element from buffer: metadata and data

As new data is provided and transferred as soon as data is read from the Queue, the get() method must not be called too frequently

class mimocorb.mimo_buffer.**Reader**(*setup_dict*)

Class to read elements from a ringbuffer (multiple-out part).

Ringbuffer elements are structured NumPy arrays and strictly **read-only**. The returned array won't change until the next Reader.get() call is performed, blocking the ringbuffer element for the time being. A program design processing the ringbuffer content has to call the Reader.get()-method in a way that minimizes the ringbuffer lock time.

methods:

- get()
- get_metadata():

data_available()

Method to check for new data and avoid blocking of consumers

get()

Get a new element from the ringbuffer. The last element obtained by calling this function is marked as "processing is done". No memory views of old elements may be accessed after calling this function (memory might change, be corrupted or be inconsistent). This function blocks if there are no new elements in the ringbuffer.

Raises

SystemExit – When the shutdown()-method of the NewBuffer object has been called, a SystemExit is raised which terminates the process.

Returns

One element (structured numpy.ndarray) of the ringbuffer as specified in the NewBuffer()-dtype-object.

Return type

numpy.ndarray

get_metadata()

Get the metadata defined for a ringbuffer element of the Reader.get()-method.

Returns

Currently a 3-tuple is returned with (counter, timestamp, deadline) which is assigned to the latest element of the ringbuffer. The content of these variables is filled by the Writer-process. The current convention is:

- counter (int): a unique, 0 based, consecutive integer referencing this element
- timestamp (float): the UTC timestamp
- **deadline (float): In a live-data environment, the dead time of the first**
writer in the analyses chain. This is meant to be the fraction of dead time to active data

capturing time (so 0.0 = no dead time whatsoever; 0.99 = only 1% of the time between this and the last element was spent with active data capturing)

Return type

tuple

class mimocorb.mimo_buffer.**Writer**(*setup_dict*)

Class to write elements into a ringbuffer (multiple-in part).

Ringbuffer elements are structured NumPy arrays. Writing is triggered by a call of `Writer.process_buffer()` or at the next call of `Writer.get_new_buffer()`. The ringbuffer element is blocked while writes to the NumPy array are permitted. A program design processing the ringbuffer content has to call the `Writer.process_buffer()` or `Writer.get_new_buffer()`-methods in a way that minimizes the ringbuffer lock time.

methods:

- `get_new_buffer()`
- `set_metadata()`
- `process_buffer()`

get_new_buffer()

Get a new free element in the ringbuffer.

The last element obtained by calling this function is marked as “ready to be processed”. No memory views of old elements may be accessed after calling this function. This function blocks if there are no free elements in the ringbuffer and always returns a valid NumPy array that can be written to.

Raises

SystemExit – When the `shutdown()`-method of the `NewBuffer` object has been called, a `SystemExit` is raised which terminates the process.

Returns

One free ringbuffer element (structured `numpy.ndarray`) as specified in the `NewBuffer()`-dtype-object. Free elements may contain older data, but they can be safely overwritten.

Return type

`numpy.ndarray`

process_buffer()

Mark the current ringbuffer element as “ready to be processed”.

The content of the array **MUST NOT** be changed after calling this function. If there is no current element, nothing happens. As the ringbuffer element is blocked while writing to the NumPy array it is recommended to call `Writer.process_buffer()` as soon as possible to minimize the ringbuffer lock time.

set_metadata(*counter, timestamp, deadtime*)

Set the metadata defined for the current ringbuffer element. If there is no current ringbuffer element (e.g. because `process_buffer()` has been called or `get_new_buffer()` has not been called yet), nothing happens. Copying metadata from a `Reader` to a `Writer` object (here called `source` and `sink`) can be done with:

```
sink.set_metadata(*source.get_metadata())
```

Parameters

- **counter** (*integer* (`np.longlong`)) – a unique, 0 based, consecutive integer referencing this element

- **timestamp** (*float (np.float64)*) – the UTC timestamp
- **deadtime** (*float (np.float64)*) – In a live-data environment, the dead time of the first writer in the analyses chain. This is meant to be the fraction of dead time to active data capturing time (so 0.0 = no dead time whatsoever; 0.99 = only 1% of the time between this and the last element was spent with active data capturing)

Collection of classes to set-up, manage and access ringbuffers and associated functions

class mimocorb.buffer_control.**buffer_control**(*buffers_dict, functions_dict, output_directory*)

Set-up and management ringbuffers and associated sub-processes

Class methods:

- **setup_buffers()**
- **setup_workers()**
- **start_workers()**
- **pause()**
- **resume()**
- **shutdown()**

display_functions()

Print list of functions and buffer associations

display_layout()

Print list of buffers

pause()

Pause data acquisition

resume()

Re-enable data acquisition after pause

setup_workers()

Set up all the (parallel) worker functions

shutdown()

Delete buffers, stop processes by calling the shutdown()-Method of the buffer manager

start_workers()

start all of the (parallel) worker functions

stop()

stop writing and reading data, allow processes to finish

class mimocorb.buffer_control.**rbExport**(*source_list=None, config_dict=None, **rb_info*)

Read data from buffer and send to requesting client (via Python yield()). Data are provided by a generator function yielding data and metadata in the `__call__()` method of the class.

class mimocorb.buffer_control.**rbImport**(*sink_list=None, config_dict=None, ufunc=None, **rb_info*)

Read data from external source (e.g. front-end device, file, simulation, etc.) and put data in mimo_buffer. Data is read by calling a user-supplied generator function for data and metadata.

class mimocorb.buffer_control.**rbObserver**(*observe_list=None, config_dict=None, **rb_info*)

Deliver data from buffer to an observer process. A tuple (data, metadata) is provided by a generator function (i.e. via yield()) implemented in the `__call__()` method of the class.

```
class mimocorb.buffer_control.rbTransfer(source_list=None, sink_list=None, config_dict=None,
                                         ufunc=None, **rb_info)
```

Read data from input buffer, filter data and write to output buffer(s) Data is provided as the argument to a user-defined filter function returning None if data is to be rejected, a number if data is to be copied to another buffer, or a list of processed input data write to additional buffers.

Args:

- buffer configurations (only one source and several sinks, no observers!)
- function ufunc() must return
 - None if data to be rejected,
 - int if only raw data to be copied to sink[0]
 - list of parameterized data to be copied to sinks[]

Action:

store accepted data in buffers

```
class mimocorb.buffer_control.rbWSObserver(observe_list=None, config_dict=None, **rb_info)
```

Deliver data from buffer to an observer process

```
class mimocorb.buffer_control.rb_toParquetfile(source_list=None, config_dict=None, **rb_info)
```

Save data a set of parquet-files packed as a tar archive

```
class mimocorb.buffer_control.rb_toTxtfile(source_list=None, config_dict=None, **rb_info)
```

Save data to file in csv-format

```
class mimocorb.buffer_control.run_mimoDAQ(setup_filename, verbose=2, debug=False)
```

Setup and run Data Acquisition suite with mimoCoRB buffer manager

The layout of ringbuffers and associated functions is defined in a configuration file in yaml format.

Functions:

- setup()
- run()
- end()

Data acquisition stops when either of the following conditions is met:

- number of requested events processed
- requested run-time reached
- input source exhausted
- end command issued from Keyboard or graphical interface

```
end(twait=3.0)
```

clean shutdown of daq suite

Arg:

twait: waiting time for processes to finish before shutdown

```
keyboard_input(cmd_queue)
```

Read keyboard input and send to Queue, running as background-thread to avoid blocking

class tc

define terminal color codes

bufferInfoGUI Graphical display of buffer status

code adapted from <https://github.com/GuenterQuast/picoDAQ>

```
mimocorb.bufferinfoGUI.bufferinfoGUI(Qcmd, Qlog, Qinfo, RBnames=['RB_1'], maxRate=100.0,  
interval=1000.0)
```

Show Buffer Manager logging messages and rate history and command buttons

Parameters

- **Qcmd** – multiprocessing.Queue for command passing to calling process
- **Qlog** – multiprocessing.Queue() for logging-info
- **Qinfo** – multiprocessing.Queue() for status info
- **RBnames** – list, buffer names, used as line labels
- **maxrate** – maximum rate for y-axis
- **interval** – update interval for graphics in ms

```
class mimocorb.bufferinfoGUI.plot_bufferinfo(Q, RBnames, maxRate=1500.0, interval=1000.0)
```

display statistics from Buffer Manager

uses multiprocessing.Queue() to display buffer information: total number of events, data acquisition rate, buffer filling level

Parameters

- **Q** – multiprocessing.Queue() for status info
- **RBnames** – list, buffer names, used as line labels
- **maxrate** – maximum rate for y-axis
- **interval** – graphics update interval in ms

plot_buffer Collection of classes with graphics functions to plot buffer data

The class `animWaveFormPlotter` is used by the class `plotBuffer()`. The `__call__()` method of this latter class is the entry point to the package.

```
class mimocorb.plot_buffer.animWaveformPlotter(conf_dict=None, dtypes=None, fig=None)
```

Oscilloscope-like display of wave from buffer data

The `__call__` method of this class receives input data and updates and redraws only the new elements of the figure created in `__init__`

init()

plot initial line objects to be animated

```
class mimocorb.plot_buffer.plot_buffer(source_list=None, sink_list=None, observe_list=None,  
config_dict=None, **rb_info)
```

Plot data using a mimiCoRB Observer

histogram_buffer collection of classes to produce histograms

Show animated histogram(s) of scalar buffer variable(s)

Because this process runs as a ‘Reader’ process, the plotting function is executed as a background task in order to avoid blocking of the main task.

The entry point is either the `__call__()` function of class `histogram_buffer`, which connects to a `mimocorb` buffer via the `rbExport` class, or directly the function `plot_Histograms()`, which receives data to be histogrammed via a multiprocessing `Queue()`.

code adapted from <https://github.com/GuenterQuast/picoDAQ>

```
class mimocorb.histogram_buffer.animHists(Hdescr, name='Histograms', fig=None)
```

display histograms, as normalised frequency distributions

```
class mimocorb.histogram_buffer.histogram_buffer(source_list=None, sink_list=None,
                                                observe_list=None, config_dict=None, **rb_info)
```

Produce Histograms of (scalar) variables in buffer.

Read data from mimiCoRB buffer using the interface class `mimo_control.rbExport` and show histograms of scalar variables selected in the configuration dictionary

Plotting is done by means of the class `plot_Histograms()` running as background process

Parameters

- **input** – configuration dictionaries
- **config_dict** – must contain a block with name ‘histograms’, formatted as
`<name>: [<min>, <max>, <nbins>, <ymin>, <ymax>, <label>, <0/1 for lin/log>]`

```
mimocorb.histogram_buffer.plot_Histograms(Q, Hdescripts, interval, name='Histograms')
```

show animated histogram(s)

Args:

- `Q`: multiprocessing.Queue()
- `Hdescripts`: list of histogram descriptors, where each descriptor is a list:
 - `min`: minimum value
 - `max`: maximum value
 - `nbins`: nubmer of bins
 - `ymax`: scale factor for highest bin ($1. = 1/Nbins$)
 - `name`: name of the quantity being histogrammed
 - `type`: 0 linear, 1 for logarithmic y scale
- `interval`: time (in s) between updates
- `name`: name of histogram window

rb_unittest: application example for `mimo_buffer`

This code may serve as a very basic starting point for own projects

Set-up: 2 ring buffers are defined:

- input Buffer `RB_1`: 10 ch x 1024 slots (int32)
- output Buffer `RB_2`: 10 ch x 2 slots/ch (float64)

simul_source: a simple template for a `mimoCoRB` source to enter simulated wave form data in a `mimoCoRB` buffer.

Input data is provided as numpy-array of shape (number_of_channels, number_of_samples).

```
simul_source.simul_source(source_list=None, sink_list=None, observe_list=None, config_dict=None,  
                           **rb_info)
```

Generate simulated data and pass data to buffer

The class `mimocorb.buffer_control.rbImport` is used to interface to the `newBuffer` and `Writer` classes of the package `mimoCoRB.mimo_buffer`

This example may serve as a template for other data sources

Parameters

config_dict – configuration dictionary

- **events_required**: number of events to be simulated or 0 for infinite
- **sleeptime**: (mean) time between events
- **random**: random time between events according to a Poission process
- **number_of_samples**, **sample_time_ns**, **pretrigger_samples** and **analogue_offset** describe the waveform data to be generated (as for oscilloscope setup)

Internal parameters of the simulated physics process (the decay of a muon) are (presently) not exposed to user.

Module **lifetime_filter**

This (rather complex) module filters waveform data to search for valid signal pulses in the channel data. The goal is to clearly identify coincidences of signals in different layers (indiating the passage of a cosmic ray particle, a muon) and find double-pulse signatures that a muon was stopped in or near a detection layer where the resulting decay-electron produced a delayed pulse. The time difference between the initial and the delayed pulses is the individual lifetime of the muon.

Wave forms passing this filter-criterion an passed on to a new buffer; the decay time and the properties of the signal pulses (height, integral and postition in time) are written to another buffer.

The relevant configuration parameters can be found in the section *calculate_decay_time*: of the configuration file.

```
lifetime_filter.calculate_decay_time(source_list=None, sink_list=None, observe_list=None,  
                                     config_dict=None, **rb_info)
```

Calculate decay time as time between double pulses

Input:

pulse wave forms

Returns:

None if failed, int or list of pulse parameters if successful

Note: output produced when filter is passed depends on number of defined sinks:

- one sink: input data
- two sinks: input data and double-pulse parameters
- three sinks: input data and double-pulse parameters separately for upwards and for downwards going decay electrons

plot: plotting waveforms from buffer using `mimoCoRB.buffer_control.ObersserverData`

```
plot_waveform.plot_waveform(source_list=None, sink_list=None, observe_list=None, config_dict=None,  
                             **rb_info)
```

Plot waveform data from mimiCoRB buffer

Parameters

input – configuration dictionary

- `plot_title`: graphics title to be shown on graph
- `min_sleeptime`: time between updates
- `sample_time_ns`, `channel_range`, `pretrigger_samples` and `analogue_offset` describe the waveform data as for oscilloscope setup

plot_histograms: histogram variable(s) from buffer using `mimoCoRB.histogram_buffer`

`plot_histograms.plot_histograms`(*source_list=None, sink_list=None, observe_list=None, config_dict=None, **rb_info*)

Online display of histogram(s) of variable(s) from mimiCoRB buffer

Parameters

input – configuration dictionary

PYTHON MODULE INDEX

I

`lifetime_filter`, 28

M

`mimocorb`, 19

`mimocorb.buffer_control`, 24

`mimocorb.bufferinfoGUI`, 26

`mimocorb.histogram_buffer`, 26

`mimocorb.mimo_buffer`, 19

`mimocorb.plot_buffer`, 26

P

`plot_histograms`, 29

`plot_waveform`, 28

R

`rb_unittest`, 27

S

`simul_source`, 27

A

`animHists` (class in `mimocorb.histogram_buffer`), 27
`animWaveformPlotter` (class in `mimocorb.plot_buffer`), 26

B

`buffer_control` (class in `mimocorb.buffer_control`), 24
`buffer_status()` (`mimocorb.mimo_buffer.NewBuffer` method), 20
`bufferinfoGUI()` (in module `mimocorb.bufferinfoGUI`), 26

C

`calculate_decay_time()` (in module `lifetime_filter`), 28

D

`data_available()` (`mimocorb.mimo_buffer.Reader` method), 22
`display_functions()` (`mimocorb.buffer_control.buffer_control` method), 24
`display_layout()` (`mimocorb.buffer_control.buffer_control` method), 24

E

`end()` (`mimocorb.buffer_control.run_mimoDAQ` method), 25

G

`get()` (`mimocorb.mimo_buffer.Observer` method), 22
`get()` (`mimocorb.mimo_buffer.Reader` method), 22
`get_metadata()` (`mimocorb.mimo_buffer.Reader` method), 22
`get_new_buffer()` (`mimocorb.mimo_buffer.Writer` method), 23

H

`histogram_buffer` (class in `mimocorb.histogram_buffer`), 27

I

`init()` (`mimocorb.plot_buffer.animWaveformPlotter` method), 26

K

`keyboard_input()` (`mimocorb.buffer_control.run_mimoDAQ` method), 25

L

`lifetime_filter` module, 28

M

`mimocorb` module, 19
`mimocorb.buffer_control` module, 24
`mimocorb.bufferinfoGUI` module, 26
`mimocorb.histogram_buffer` module, 26
`mimocorb.mimo_buffer` module, 19
`mimocorb.plot_buffer` module, 26
module
 `lifetime_filter`, 28
 `mimocorb`, 19
 `mimocorb.buffer_control`, 24
 `mimocorb.bufferinfoGUI`, 26
 `mimocorb.histogram_buffer`, 26
 `mimocorb.mimo_buffer`, 19
 `mimocorb.plot_buffer`, 26
 `plot_histograms`, 29
 `plot_waveform`, 28
 `rb_unittest`, 27
 `simul_source`, 27

N

`new_observer()` (`mimocorb.mimo_buffer.NewBuffer` method), 20

`new_reader_group()` (*mimocorb.mimo_buffer.NewBuffer* method), 21

`new_writer()` (*mimocorb.mimo_buffer.NewBuffer* method), 21

`NewBuffer` (class in *mimocorb.mimo_buffer*), 19

O

`Observer` (class in *mimocorb.mimo_buffer*), 21

P

`pause()` (*mimocorb.buffer_control.buffer_control* method), 24

`pause()` (*mimocorb.mimo_buffer.NewBuffer* method), 21

`plot_buffer` (class in *mimocorb.plot_buffer*), 26

`plot_bufferinfo` (class in *mimocorb.bufferinfoGUI*), 26

`plot_histograms` module, 29

`plot_Histograms()` (in module *mimocorb.histogram_buffer*), 27

`plot_histograms()` (in module *plot_histograms*), 29

`plot_waveform` module, 28

`plot_waveform()` (in module *plot_waveform*), 28

`process_buffer()` (*mimocorb.mimo_buffer.Writer* method), 23

R

`rb_toParquetfile` (class in *mimocorb.buffer_control*), 25

`rb_toTxtfile` (class in *mimocorb.buffer_control*), 25

`rb_unittest` module, 27

`rbExport` (class in *mimocorb.buffer_control*), 24

`rbImport` (class in *mimocorb.buffer_control*), 24

`rbObserver` (class in *mimocorb.buffer_control*), 24

`rbTransfer` (class in *mimocorb.buffer_control*), 24

`rbWSObserver` (class in *mimocorb.buffer_control*), 25

`Reader` (class in *mimocorb.mimo_buffer*), 22

`resume()` (*mimocorb.buffer_control.buffer_control* method), 24

`resume()` (*mimocorb.mimo_buffer.NewBuffer* method), 21

`run_mimoDAQ` (class in *mimocorb.buffer_control*), 25

`run_mimoDAQ.tc` (class in *mimocorb.buffer_control*), 25

S

`set_ending()` (*mimocorb.mimo_buffer.NewBuffer* method), 21

`set_metadata()` (*mimocorb.mimo_buffer.Writer* method), 23

`setup_workers()` (*mimocorb.buffer_control.buffer_control* method), 24

`shutdown()` (*mimocorb.buffer_control.buffer_control* method), 24

`shutdown()` (*mimocorb.mimo_buffer.NewBuffer* method), 21

`simul_source` module, 27

`simul_source()` (in module *simul_source*), 27

`start_workers()` (*mimocorb.buffer_control.buffer_control* method), 24

`stop()` (*mimocorb.buffer_control.buffer_control* method), 24

W

`Writer` (class in *mimocorb.mimo_buffer*), 23