

---

**mimoCoRB**

***Release 1.0.2rc1***

**C. Mayer, K. Heitlinger, G. Quast**

**Sep 02, 2024**



**CONTENTS:**

- 1 mimoCoRB - multiple-in multile-out Configurable Ring Buffer: Overview 1**
  - 1.1 For Developers: Description of components . . . . . 1
  - 1.2 Simple application example . . . . . 2
- 2 For application developers: Access Classes in the module *buffer\_control* 7**
  - 2.1 Application examples . . . . . 10
- 3 Module Documentation 21**
- Python Module Index 33**
- Index 35**



---

## MIMOCORB - MULTIPLE-IN MULTIPLE-OUT CONFIGURABLE RING BUFFER: OVERVIEW

**mimoCoRB**: multiple-in multiple-out Configurable Ring Buffer

The package **mimoCoRB** provides a central component of each data acquisition system needed to record and pre-analyse data from randomly occurring processes. Typical examples are waveform data as provided by detectors common in quantum mechanical measurements, or in nuclear physics, particle physics and astro particle physics, e. g. photo tubes, Geiger counters, avalanche photo-diodes or modern SiPMs.

The random nature of such processes and the need to keep read-out dead times low requires an input buffer for fast collection of data and an efficient buffer manager delivering a constant data stream to the subsequent processing steps. While a data source feeds data into the buffer, consumer processes receive the data to filter, reduce, analyze or simply visualize the recorded data. In order to optimally use the available resources, multi-core and multi-processing techniques must be applied. Data consumers may be obligatory ones, i. e. data acquisition pauses if all input buffers are full and an obligatory consumer is still busy processing. A second type of consumers (random consumers or “observers”) receive an event copy from the buffer manager upon request, without interrupting the data acquisition process. Typical examples of random consumers are displays of a subset of the waveforms or of intermediate analysis results.

This project originated from an effort to structure and generalize data acquisition for several experiments in advanced physics laboratory courses at Karlsruhe Institute of Technology (KIT) and has been extensively tested with Ubuntu Linux.

As a simple, stand-alone demonstration, we provide simulated signals as would be produced by a detector for cosmic muons with four detection layers. Occasionally, such muons stop in an absorber between the 2nd and 3rd layer, where they decay at rest and emit a high-energetic electron recorded as a 2nd pulse in one or two of the detection layers. After data acquisition, a search for typical pulse shapes is performed and data with detected double pulses are selected and copied into a second buffer. A third buffer receives data in a reduced format which only contains the parameters of accepted pulses. These data and the waveforms of all double-pulses are finally stored on disk. Such an application is a very typical example of the general process of on-line data processing in modern physics experiments and may serve as a starting point for own projects.

### 1.1 For Developers: Description of components

The following paragraphs provide some insight into the inner working of the components of *mimoCoRB* for interested users and for developers wanting to help improving the package. Application developers should use the more convenient access classes described below to build an application based on the *mimoCoRB* framework.

In order to decouple the random occurrence of “events” at the start of a data-acquisition chain one needs a buffer capable of rapidly storing new incoming data and delivering a constant data stream to subsequent consumer processes. This is typically implemented as a first-in, first out ringbuffer providing storage space in memory for incoming data, which is released and overwritten by new data when all consuming processes have finished.

As digital filtering of incoming data may be very CPU intensive, multi-processing and multi-core capable components are needed to ensure sufficient compute power to process and analyze all data. *mimoCoRB.mimo\_buffer* implements such a buffer allowing multiple processes to read (“multiple out”) or write (“multiple in”) to a shared buffer space. Access to common memory storage and the synchronization of processes is achieved using the *shared\_memory*, *Process* and *Queue* modules from the *multiprocessing* package.

Because processing of the data, i. e. digital filtering, selection, compression and storage or real-time visualization of the data can be a complex workflow, data buffers may be arranged in chains where one or several reader processes associated to a buffer write to one or several output buffer(s).

The central component takes care of memory management and access control provided by the class **newBuffer**. To control the data flow in a full data acquisition suite, three types of access are foreseen, implemented as **Writer**, **Reader** and **Observer** classes. Readers of the same type are grouped together for multi-processing of compute-intense tasks and form a Reader-group. Observers receive only a sub-set of the data and are mainly intended to be used for visual inspection or graphical representation of samples of the recorded or processed data.

Processes for data provisioning from front-end hardware or from other sources, like disk files, web streams or simulation, rely on the Writer class; similarly, processes reading data for filtering and transfer to subsequent buffers or to analyse, transform or extract data to storage media use the Reader class. Note that the buffer manager ensures that every slot assigned to a Reader (or a group of Readers) is actually processed; therefore, input to a buffer blocks if the buffer is filled up completely. The Writer class resumes data input as soon as a Reader or member of a Reader-group has finished processing and thus freed a slot in the buffer.

Multiprocessing is enabled by use of the *shared\_memory* module of the *multiprocessing* package available since Python 3.8 for direct access to shared memory across processes. Other modules of the package (*Process*, *Lock*, *Event*, and *SimpleQueue* or *Queue*) are used to create and control sub-processes and for signaling and message or data exchange across processes.

The format of data stored in the buffers is based on structured *numpy* arrays with (configurable) field names and *numpy* *dtypes*. Each buffer entry is also associated with a unique event number, a time stamp and a deadtime fraction to be provided by the initial data producer. The deadtime accounts for inefficiencies of the data acquisition due to processing in *mimoCoRB*. These metadata are set by the initial producer and must not be changed at a later stage in the processing chain.

## 1.2 Simple application example

An application example of *mimo\_buffer* is shown below; it is also provided as a unit test. The set-up is as follows:

Two ring buffers are defined:

- input Buffer RB\_1: 10 ch x 1024 slots (int32)
- output Buffer RB\_2: 10 ch x 2 slots/ch (float64)

Simple data is filled into RB\_1, copied and extended by a process writing data into RB\_2, and finally a reader process to check integrity and completeness of the data. The most complex part of the code is in function *run\_control()*, which demonstrates how to set up the buffers, define Reader and Writer instances and start the parallel processes for generating, processing and reading the data.

The example including comment lines for explanation is shown here:

```
import time
import unittest
import numpy as np
from multiprocessing import Process, Value
from mimocorb import mimo_buffer as bm
```

(continues on next page)

(continued from previous page)

```

# global variables
N_requested = 1000 # number of data injections ("events")
Time_tick = 0.001 # time between events
Ncpul = 2 # number of parallel analyzer processes

def data_generator(sink_dict):
    """writes continuously rising integers to buffer specified in sink_dict
    """
    sink = bm.Writer(sink_dict)
    n=0
    # inject data
    for x in range(N_requested):
        buffer = sink.get_new_buffer() # get new buffer and pass last item
        # random wait for next data item
        time.sleep(-Time_tick*np.log(np.random.rand() ))
        # fill "data"
        n += 1
        buffer[:] = n
    # process last data item
    sink.process_buffer()

def analyzer(source_dict, sink_dict):
    """read from source and write first element and a time difference to sink
    """
    source = bm.Reader(source_dict)
    sink = bm.Writer(sink_dict)
    start_time = time.time()

    while True:
        input_data = source.get()
        output_data = sink.get_new_buffer()
        # process data
        output_data[0] = input_data[0]
        # mimick processing time
        time.sleep(2*Time_tick)
        output_data[1] = time.time() - start_time

        #
        sink.process_buffer()

def check_result(source_dict, res):
    """reads RB_2 and sum up the integer content

    sum is returned as shared memory Value-object
    """
    source = bm.Reader(source_dict)
    sum_rb = 0
    while True:
        input_data = source.get()
        res.value +=int(input_data[0])

def run_control():

```

(continues on next page)

(continued from previous page)

```

"""Setup buffers, start processes and shut_down when 1st writer done
"""

# Create ring buffers: #2: 10 channel, 2 value per channel
# (1: buffer content; 2: time difference as int)
# d_type = [('chA', np.float)] #not necessary: always the same type
generator_buffer = bm.NewBuffer(10, 1, np.int32)
eval_buffer = bm.NewBuffer(10, 2, np.float32)

# create readers first
source_dic_gen = generator_buffer.new_reader_group()
source_dic_eval = eval_buffer.new_reader_group()

# Create worker processes (correct sequence: first action as last)
process_list = []
# evaluation to test ring buffer behavior
result = Value('i', 0) # int variable in shared meomry
process_list.append(Process(target=check_result,
                           args=(source_dic_eval, result)))
# data transfer between the 2 buffers: generator_buffer -> eval_buffer
sink_dic_eval = eval_buffer.new_writer()
# work with all cpu's requested
number_of_workers = Ncpu1
for i in range(number_of_workers):
    process_list.append(Process(target=analyzer,
                              args=(source_dic_gen, sink_dic_eval)))

# fill buffer (generator_buffer) with data first
sink_dic_gen = generator_buffer.new_writer()
process_list.append(Process(target=data_generator,
                           args=(sink_dic_gen,)))

for p in process_list:
    p.start()

run_active = True
while run_active:
    run_active = False if process_list[-1].exitcode==0 else True
    time.sleep(0.1) # wait
time.sleep(0.1) # some grace-time for readers to finish

generator_buffer.shutdown()
eval_buffer.shutdown()
del generator_buffer, eval_buffer

for p in process_list:
    p.join()

return result.value

class RPTest(unittest.TestCase):

```

(continues on next page)



(continued from previous page)

```
def test_process(self):
    # start python test module and check result
    a = run_control()
    expected_result = N_requested*(N_requested+1)//2
    # expected result: sum(i); i = 1, N_requested
    self.assertEqual(a, expected_result)

if __name__ == "__main__":
    unittest.main(verbosity=2)
```



---

## FOR APPLICATION DEVELOPERS: ACCESS CLASSES IN THE MODULE *BUFFER\_CONTROL*

To facilitate user interaction with the buffer manager, a set of additional classes is provided in the module *buffer\_control* to set-up and manage cascades of ringbuffers and the associated functions for filling, filtering and extracting data. These classes are most interesting for application developers wanting to build upon the *mimoCoRB* framework.

The classes are:

- **class *buffer\_control*:**  
Set-up and management of ringbuffers and associated sub-processes. This is the overarching class with access to all created buffers and sub-processes.
- **class *rbImport*:**  
Read data from source (e.g. a front-end like a PicoScope USB oscilloscope, or from a file or simulation) and import data and metadata in a *mimo\_buffer* by calling user-supplied Python generator (i.e. via *'yield'*). In this approach, *mimiCoRB* “pulls” data.
- **class *rbPut*:**  
Read data from source (e.g. a front-end like a PicoScope USB oscilloscope, or from a file or simulation) and put data in a *mimo\_buffer* by calling a Python function, thus pushing the data under control of the reading application. This method is useful in cases where the application providing the data has its own event loop (driving e. g. its own graphical user interface).
- **class *rbTransfer*:**  
Read data from a *mimo\_buffer*, filter and/or reformat data and write to output *mimo\_buffer(s)*. Data is provided as the argument to a user-defined filter function, returning *None* if data is to be discarded, a number if data is to be copied to another buffer, or - optionally - a list of transformed data records produced from processed input data. If such data are provided, a respective number of ringbuffers as destination must be configured.
- **class *rbExport*:**  
Read data from *mimo\_buffer* and analyze (with user-supplied code), without writing to another ringbuffer. Data are provided to the user function as a Python generator in the *\_\_call\_\_()* method of the class yielding a tuple of data and metadata.
- **class *rbObserver***  
Deliver data from a buffer to an observer process. A tuple (data, metadata) is provided by a Python generator implemented in the *\_\_call\_\_()* method of the class.
- **class *rb\_toTxtfile*:**  
Save *mimo\_buffer* data to a file in csv-format. The header line of this file contains the keys of the respective columns, which are derived from the datatype of the structured ringbuffer array. Aliases for improved clarity can be provided in the configuration file.
- **class *rb\_toParquetfile*:**  
Save *mimo\_buffer* data to an archive in *tar* format; each data record is packed in *Parquet* format.

- **class *run\_mimoDAQ*:**


Setup and run a Data Acquisition suite with the mimoCoRB buffer manager. The layout of ringbuffers and associated functions are defined in a configuration file in *yaml* format. All configured functions are executed as worker processes in separate sub-processes and therefore optimal use is made of multi-core architectures.

- **class *bufferinfoGUI*:**

**A graphical interface showing buffer rates and status information and**  
providing some control buttons interacting with the *run\_mimoDAQ* class.


These classes shield much of the internal complexity from the user, who can thus concentrate on writing the pieces of code needed to acquire and process the data. The access classes expect as input lists of dictionaries with the parameters of buffers to read from (**source\_list**), to write to (**sink\_list**) or to observe (**observe\_list**). An additional dictionary (**config\_dict**) provides the parameters needed for the specific functionality, for example names of functions to read, filter or manipulate data or the names of target files. The interface for passing data between the user-defined functions and ringbuffers relies on Python generators (i.e. the *yield* instruction).

The overarching class **buffer\_control** provides methods to setup buffers and worker processes and to control the data acquisition process. The methods collected in the class *run\_mimoDAQ*, in particular the function **run\_mimoDAQ**, contains the code needed to run an instance of a data-acquisition suite defined in a configuration file specifying the ring buffers and associated, user-defined functions for data provisioning, filtering and storage. *run\_mimoDAQ* is controlled either by keyboard commands or from a graphical user interface; pre-defined conditions on the total number of events processed, the duration of the data taking run or finishing of the writer process to the first buffer due to source exhaustion can also be defined to end data taking. The class structure and dependencies are shown in the figure below.

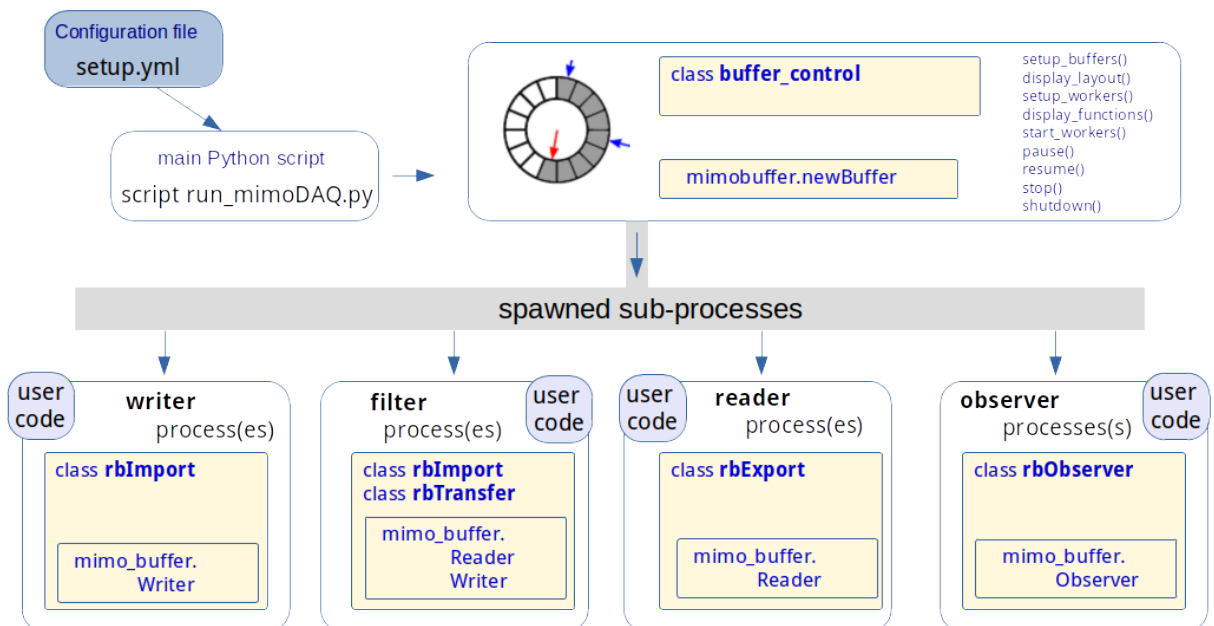

A sequence diagram of a typical application, shown below, illustrates the interplay and dependencies of the classes described above. The script *run\_daq.py* creates an instance of *run\_mimoDAQ* and starts its *run()*-method. The interaction with the user-supplied filter functions is handled by methods of the class *buffer\_control*.

## mimocorb sequence (main parts)

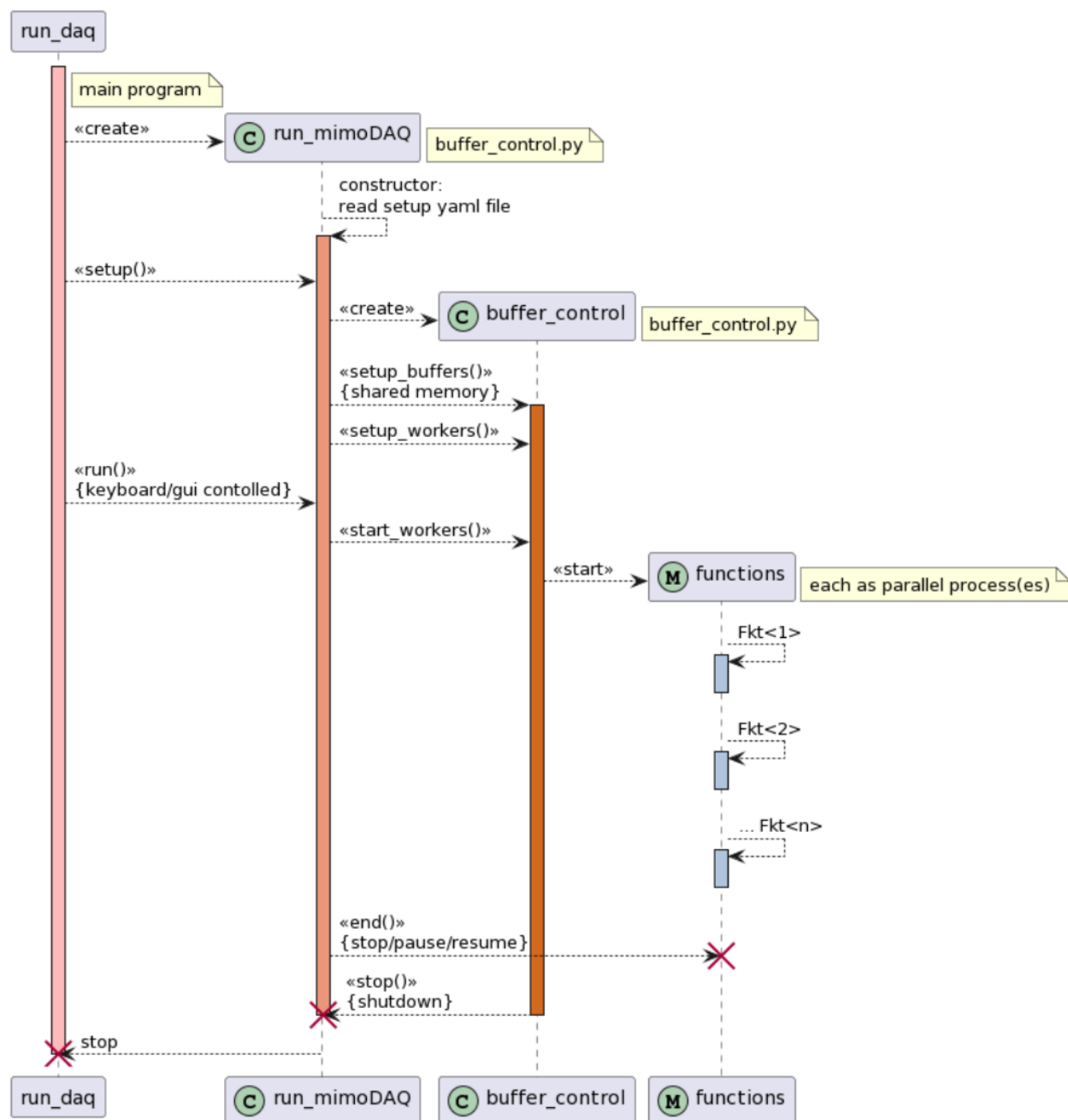

For complex setups and longer data-taking periods it is important to gain a quick overview of the status of all buffers and to monitor long-term stability. Therefore, a graphical display with the processing rate of all buffers is provided by the class **bufferinfoGUI**. A text window receives frequent updates of the number of events processed by each buffer and of the buffer fill-levels. Clickable control buttons send information via a dedicated command queue to the calling process *run\_mimoDAQ* and enable pausing, resuming and controlled ending of the data-acquisition processes.

The suggested structure of the project work-space for *mimoCoRB* applications is as follows:

```
|--> <user working directory>          # the main configuration script resides here
```

(continues on next page)

(continued from previous page)

```
| --> modules      # project-specific, user-supplied python code
| --> config        # configuration files in yaml format
| --> target        # output of data-acquisition run(s)
```

For illustration and as a starting point for own applications, a stand-alone example is provided as part of the package, as described in the following section.

## 2.1 Application examples

The subdirectory *examples/* contains some complete application use cases based on input waveforms that are generated by a simulator in real-time.

The central piece of every *mimoCoRB* application is the configuration file; examples of different complexity are provided in *examples/\*\_setup.yaml*. Code snippets for data input, filtering and output as well as configuration files are provided in the subdirectories *examples/modules/* and *examples/config/*, respectively.

Waveform data, as provided by, for example, a multi-channel digital oscilloscope, are generated and filled into the first one of a cascaded set of ringbuffers and passed on to subsequent buffer stages depending on filter conditions. Data in the last buffer are recorded to disk. The configuration files and the recorded data files are stored in the subdirectory *examples/target/<projectname>\_<date\_and\_time>*.

All examples run stand-alone and use as input simulated waveform data of short pulses as they arise e. g. in scintillation or semi-conductor detection devices for single particles. The simulated physics process corresponds to signatures produced by cosmic muons penetrating several layers of scintillator material.

### Simple Example

A very simple example consists of recording two input channels from two redundant detectors. The buffer configuration is defined in the file *examples/demo\_setup.yaml*, shown here:

```
# Configuration for recording two channels with mimoCoRB
# -----
#
# illustrates the structure of a mimoCoRB setup file:
# the blocks
# - RingBuffer:
# - Functions:
# must be present.
# The names RB_i and FKT_i with i = 1, 2, 3, ... are fixed.
# Optionally, a block
# - FunctionConfigs:
# can be specified to provide the configuration for each one of the
# assigned functions. These may also be specified in a separate
# configuration file for all functions, or in dedicated files for
# each function using the keys "config_file:" in the FTK_main section
# or in the FKT_i: sections, respectively.
#
# In this example, two buffers are configured:
# - RB_1 for raw waveforms
# - RB_2 for derived pulse parameters
#
# Data from RB_2, the result buffer, are saved to a file in csv (text) format.
# Data from RB_1 are also passed to an observer process driving a real-time display,
```

(continues on next page)

(continued from previous page)

```
# and are also passed to a Reader process driving a real-time histogram display.
#
# Notes:
#
# 1. additional config files controlling the user functions are
#    located in the subdirectory config/
#
# 2. necessary application-specific user code is located
#    in the subdirectory modules/
#
# -----
#
```

**RingBuffer:**

```
# define ring buffers
- RB_1:
  # raw input data buffer (from picoScope, file or simulation)
  number_of_slots: 16
  channel_per_slot: 500
  data_type:
    1: ['chA', "float32"]
    2: ['chB', "float32"]
- RB_2:
  # buffer with correct signature double pulse parameters
  number_of_slots: 16
  channel_per_slot: 1
  data_type:
  data_type:
    1: ['chA_height', "float32"]
    2: ['chA_position', "int32"]
    3: ['chA_integral', "float32"]
    4: ['chB_height', "float32"]
    5: ['chB_position', "int32"]
    6: ['chB_integral', "float32"]
```

**Functions:**

```
# define functions and ringbuffer assignment

- Fkt_main:
  # runtime: 60 # desired runtime in seconds
  runevents: 1000
  config_file: "config/spectrum_config.yaml"

- Fkt_1:
  ## for simulation source
  file_name: "modules/simul_source"
  fkt_name: "simulation_source"
  ## for data from file
  #file_name: "modules/file_source"
  #fkt_name: "tar_parquet_source"
  num_process: 1
  RB_assign:
```

(continues on next page)

(continued from previous page)

```

    RB_1: "write"

- Fkt_2:
  file_name: "modules/exporters"
  # fkt_name: "save_parquet"
  fkt_name: "drain"
  num_process: 1
  RB_assign:
    RB_1: "read"      # waveform to save

- Fkt_3:
  file_name: "modules/spectrum_filter"
  fkt_name: "find_peaks"
  num_process: 2
  RB_assign:
    RB_1: "read"
    RB_2: "write"

- Fkt_4:
  file_name: "modules/exporters"
  # fkt_name: "save_to_txt" # save data to text
  fkt_name: "drain"        # no saving of data
  num_process: 1
  RB_assign:
    RB_2: "read"

```

The example coming with this package contains two more convenience functions, one for an observer process displaying a random sample of waveforms in an oscilloscope display, and a second one for on-line analysis and histogramming of buffer data. The necessary addendum to the configuration looks as follows:

```

# --- the following functions are optional

- Fkt_5:
  file_name: "modules/plot_waveform"
  fkt_name: "plot_waveform"
  num_process: 1
  RB_assign:
    RB_1: "observe"

- Fkt_6:
  file_name: "modules/plot_histograms"
  fkt_name: "plot_histograms"
  num_process: 1
  RB_assign:
    RB_2: "read" # pulse parameters

```

These additional functions cover very general use cases and rely on the modules *mimocorb.plot\_buffer* and *mimocorb.histogram\_buffer*, which provide animated displays of waveforms similar to an oscilloscope and a histogram package for life-updates of frequency distributions of scalar variables.

The last *yaml* block contains configuration parameters needed for some of the functions.

```

FunctionConfigs:
# configuration of functions provided by mimoCoRB

```

(continues on next page)



(continued from previous page)

```

save_to_txt:
    filename: "spectrum"

save_parquet:
    filename: "spectrum"

drain:      # has no parameters, empty field is sufficient
key:

```

As an alternative, configuration parameters for functions can also be provided in a separate *yaml* file in the block *FKT\_main*. (line *config\_file*: "*config/spectrum\_config.yaml*"). It is also possible to specify configuration files for each of the functions individually in the respective block *FKT\_i*, *i*=1,2,...

The declared functions must support a unique calling interface and are started as sub-processes. Lists of dictionaries provide the necessary information to connect to the buffer manager via the *Writer*, *Reader* or *Observer* classes of the package. This information comprises the pointer to the shared buffer manager as well as pointers to instances of the functions *Event()* or *Queue()* from the multiprocessing package to enable communication and data transfers across processes. A further dictionary (*config\_dict*) provides the function-specific configuration parameters discussed previously. The keyword dictionary *rb\_info* specifies whether writer, reader or observer functionality is required. It contains a copy of the ring-buffer assignment block („RB\_assign:") from the main setup file on function level. Its purpose is to facilitate the ring-buffer access part within a function or class.

The function interface looks as follows:

```

def <function_name>(  
    source_list=None, sink_list=None, observe_list=None, config_dict=None, **rb_info):

```

This interface must be respected by any user function. The argument list must also be passed to instances of the access classes *rbImort*, *rbExport*, *rbTransfer* or *rbObserve*. An example of a user function in the directory *modules/* to write buffer data to a text file is shown below:

```

"""Module exporters to handle file I/O for data in txt and parquet format
    This module relies on classes in mimocorb.buffer_control
"""
from mimocorb.buffer_control import rb_toTxtfile, rb_toParquetfile
def save_to_txt(source_list=None, sink_list=None, observe_list=None, config_dict=None,
    **rb_info):
    sv = rb_toTxtfile(source_list=source_list, config_dict=config_dict, **rb_info)
    sv()

```

Running the example with the command

```
../run_daq.py demo_setup.yaml
```

yields the following output on screen:

```

*==* script ../run_daq.py running

2 buffers created... List of buffers
RB_1 16 500
RB_2 16 1
List of functions
FKT_1 simul_source (1)  {'RB_1': 'write'}
FKT_2 find_peaks (2)   {'RB_1': 'read', 'RB_2': 'write'}
FKT_3 save_to_txt (1)  {'RB_2': 'read'}

```

(continues on next page)

(continued from previous page)

```
FKT_4 plot_waveform (1)    {'RB_1': 'observe'}
FKT_5 plot_histograms (1)  {'RB_2': 'read'}
```

Two buffers are created in this case, *RB1\_1* and *RB\_2*. *RB\_1* with two channels with 500 samples each is the input buffer. *RB\_2* contains 6 scalar variables and is the output buffer, which is filled by the function *find\_peaks* with two active workers. The functions *save\_to\_txt* and *plot\_histograms* read from this buffer and store data to disk or show histograms, respectively. The function *plot\_waveform* takes random samples from *RB\_1* and displays the raw waveform data.

This example serves as a convenient starting point for own application development. The code in *simul\_source.py*, shown below, is a very general example for data input. Only the function *pulseSimulator()* needs to be replaced by a function providing data from your own source. The code is shown here:

```
def simulation_source(source_list=None, sink_list=None, observe_list=None, config_
dict=None, **rb_info):
    """
    General example for data import from external source
    (here: generation of simulated data with module pulseSimulator)

    Uses class mimocorb.buffer_control/rbImport to interface to the
    newBuffer and Writer classes of the package mimoCoRB.mimo_buffer

    mimiCoRB interacts with this code via a generator (*yield_data()*),
    which itself received data via the *__call__* function of the class
    *dataSource* providing the input data. Configuration parametes
    in the dictionary *config_dict* are passed to this class during
    initialisation. Parameters of the configured buffers are set after
    after initialisation.

    This example may serve as a template for other data sources
    """

    # define and instantiate external data source
    source = dataSource(config_dict)

    def yield_data():
        """generate simulated data, called by instance of class mimoCoRB.rbImport"""
        event_count = 0
        while True:
            data = source()
            # deliver pulse data (and no metadata - these are added by rbImport)
            yield (data, None)
            event_count += 1

    # get buffer configuration
    sink_dict = sink_list[0]
    number_of_channels = len(sink_dict["dtype"])
    number_of_values = sink_dict["values_per_slot"]
    channel_names = [sink_dict["dtype"][i][0] for i in range(number_of_channels)]
    # consistency check
    if "number_of_samples" not in config_dict:
        pass
    else:
```

(continues on next page)

(continued from previous page)

```

    if number_of_values != config_dict["number_of_samples"]:
        print("! Config Error: requested number of samples does not match buffer_
↪size !")
        sys.exit("requested number of samples does not match buffer size !")
    source.init(number_of_channels, number_of_values, channel_names)

    # instantiate buffer manager interface
    rbImporter = rbImport(config_dict=config_dict, sink_list=sink_list, ufunc=yield_data,
↪ **rb_info)
    rbImporter.logger.info(f"sub-process {__name__} started")
    # start __call__ method of rbImport instance
    rbImporter()


```

### Complex Example

In the second, more complex example discussed here we consider multiple pulses on the same channels, where, e. g. the first pulse originates from a detected muon and the second, later one from a decay electron of a muon that has been stopped in or near a detection layer.

The raw data are analyzed, and accepted data with a double-pulse signature are selected and directly passed on to a second ring buffer. A third buffer contains only the information on found signal pulses; a result file in *csv* format contains the data extracted from this buffer.

A graphical representation of the set-up is shown in the figure below [source: Master's Thesis Christoph Mayer, ETP 2022]. Note that the oscilloscope as the data-acquisition front-end is replaced by a signal simulation in the example provided.

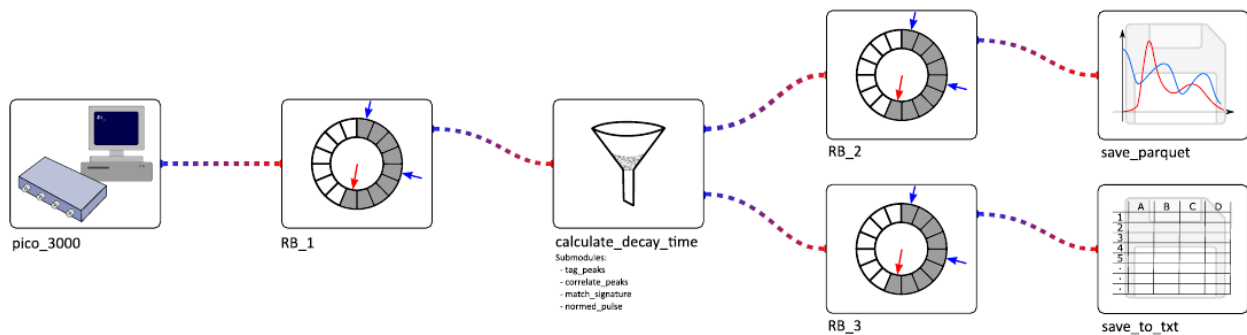

The buffer layout and the associated functions are defined in the main configuration file *lifetime\_setup.py*, which serves as the input to the execution script *run\_daq.py* in the top-level directory of the package. The *python* files *simul\_source.py*, *lifetime\_filter.py* and *exporters.py* contain the user code for data generation, analysis and filtering and extraction of the finally accepted data to disk files. The *.yaml* files *simulation\_config.yaml* and *save\_lifetimes.yaml* contain configurable parameters provided to these functions.

This example is executed from the directory *examples/* by entering:

```
../run_daq.py lifetime_setup.yaml
```

Again, the screen output gives an overview of the generated buffers and the functions writing to and reading from them:

```

*** script ../run_daq.py running

3 buffers created... List of buffers
RB_1 128 4250
RB_2 128 4250

```

(continues on next page)

(continued from previous page)

```

RB_3 32 1
List of functions
FKT_1 simul_source (1)  {'RB_1': 'write'}
FKT_2 calculate_decay_time (2)  {'RB_1': 'read', 'RB_2': 'write', 'RB_3': 'write'}
FKT_3 save_to_txt (1)  {'RB_3': 'read'}
FKT_4 save_parquet (1)  {'RB_2': 'read'}
FKT_5 plot_waveform (1)  {'RB_2': 'observe'}
FKT_6 plot_histograms (1)  {'RB_3': 'read'}

```

The input *yaml* file for this example looks as follows:

```

# Application example for mimoCoRB
# -----
#
# three buffers:
# - RB_1 for (simulated) raw waveforms
# - RB_2 for selected double-pulses
# - RB_3 for derived pulse parameters
#
# data from RB_2 and RB_3 are saved to files in tarred parquet format
# or in text format.
#
# data from RB_2 are passed to an observer process driving a real-time display
#
# Notes:
#
# 1. additional config files controlling the user functions are
#    located in the subdirectory config/
# 2. user necessary application-specific user code is located
#    in the subdirectory modules/
# -----
#
# general control options
output_directory: target      # directory to store output data
GUI_control: true             # control daq via Grapical User Interface
KBD_control: true             # control daq via KeyBoard

RingBuffer:
  # define ring buffers
  - RB_1:
      # raw input data buffer (waveforms from PicoScope, file_source or simulation)
      number_of_slots: 128
      channel_per_slot: 4250
      data_type:
        1: ['chA', "float32"]
        2: ['chB', "float32"]
        3: ['chC', "float32"]
        4: ['chD', "float32"]
  - RB_2:
      # buffer with accepted signatures (here double-pulses)

```

(continues on next page)

(continued from previous page)

```

number_of_slots: 128
channel_per_slot: 4250
data_type:
  1: ['chA', "float32"]
  2: ['chB', "float32"]
  3: ['chC', "float32"]
  4: ['chD', "float32"]
- RB_3:
  # buffer with pulse parameters (derived from waveforms)
  number_of_slots: 32
  channel_per_slot: 1
  data_type:
    1: ['decay_time', "int32"]
    3: ['1st_chA_h', "float32"]
    4: ['1st_chB_h', "float32"]
    5: ['1st_chC_h', "float32"]
    6: ['1st_chA_p', "int32"]
    7: ['1st_chB_p', "int32"]
    8: ['1st_chC_p', "int32"]
    9: ['1st_chA_int', "float32"]
    10: ['1st_chB_int', "float32"]
    11: ['1st_chC_int', "float32"]
    12: ['2nd_chA_h', "float32"]
    13: ['2nd_chB_h', "float32"]
    14: ['2nd_chC_h', "float32"]
    15: ['2nd_chA_p', "int32"]
    16: ['2nd_chB_p', "int32"]
    17: ['2nd_chC_p', "int32"]
    18: ['2nd_chA_int', "float32"]
    19: ['2nd_chB_int', "float32"]
    20: ['2nd_chC_int', "float32"]
    21: ['1st_chD_h', "float32"]
    22: ['1st_chD_p', "int32"]
    23: ['1st_chD_int', "float32"]
    24: ['2nd_chD_h', "float32"]
    25: ['2nd_chD_p', "int32"]
    26: ['2nd_chD_int', "float32"]
Functions:
  # define functions and assignments
  - Fkt_main:
    config_file: "config/simulation_config.yaml"
  - Fkt_1:
    file_name: "modules/simul_source"
    fkt_name: "simul_source"
    num_process: 1
    RB_assign:
      RB_1: "write"
  - Fkt_2:
    file_name: "modules/lifetime_filter"
    fkt_name: "calculate_decay_time"
    num_process: 2
    RB_assign:

```

(continues on next page)

(continued from previous page)

```

        RB_1: "read"      # input
        RB_2: "write"     # waveform to save (if double pulse was found)
        RB_3: "write"     # pulse data
- Fkt_3:
    file_name: "modules/exporters"
    fkt_name: "save_to_txt"
    config_file: "config/save_lifetime.yaml"
    num_process: 1
    RB_assign:
        RB_3: "read"      # pulse data
- Fkt_4:
    file_name: "modules/exporters"
    fkt_name: "save_parquet"
    num_process: 1
    RB_assign:
        RB_2: "read"      # waveform to save

```

The functions for plotting result variables and for histogramming are not shown, but also contained in the example configuration provided as part of the package, in the same way as explained for the first example.

In addition to the remarks concerning configuration files, it is worth noticing here that a dedicated configuration file is specified in a separate *yaml* file for function *Fkt\_3*. This possibility is particularly useful if the same function code is used to handle data from different buffers, e.g. the file name and special formatting for writing buffer contents to a file in *csv* format.

### Even more complex example

A similar, but even more complex case with two streams going to different output buffers and files is specified in the file *spin\_setup.yaml*. Starting this example with `> ./run_daq.py spin_setup.yaml` leads to the following terminal output:

```

*** script ../run_daq.py running

4 buffers created... List of buffers
RB_1 128 4250
RB_2 128 4250
RB_3 32 1
RB_4 32 1
List of functions
FKT_1 simul_source (1) {'RB_1': 'write'}
FKT_2 calculate_decay_time (2) {'RB_1': 'read', 'RB_2': 'write', 'RB_3': 'write', 'RB_
↪4': 'write'}
FKT_3 save_to_txt (1) {'RB_3': 'read'}
FKT_4 save_to_txt (1) {'RB_4': 'read'}
FKT_5 save_parquet (1) {'RB_2': 'read'}
FKT_6 plot_waveform (1) {'RB_2': 'observe'}
FKT_7 plot_histograms (1) {'RB_3': 'read'}
FKT_8 plot_histograms (1) {'RB_4': 'read'}

Running 260s RB_1: 76409(0) 302Hz RB_2: 6118(0) 19Hz RB_3: 2952(1) 9.99Hz RB_4: 3166(1)
↪8.99Hz

```

A screenshot of a data-acquisition run with input from simulated data is shown in the figure below.







## MODULE DOCUMENTATION

### **mimo-ringbuffer:**

Module implementing a multiple-in multiple-out ringbuffer appropriate for multi-core multiprocessing.

The ringbuffer creation and management is handled by the class `NewBuffer`. Buffer access is managed by the `Reader`, `Writer` and `Observer` classes.

classes:

- `NewBuffer`: create a new ringbuffer, assign writer(s) and reader(s) or observer(s)

methods:

- `new_writer`
  - `new_reader_group`
  - `new_observer`
  - `buffer_staus`
- `Writer`: write elements into a ringbuffer
  - `Reader`: read all elements from a ringbuffer
  - `Observer`: read selected elements from a ringbuffer.

**class** `mimocorb.mimo_buffer.NewBuffer`(*number\_of\_slots, values\_per\_slot, dtype, debug=False*)

Class to create a new ringbuffer object according to the 'FIFO' principle (first-in first-out).

Memory shares, IPC queues, lock and event objects as well as background threads are defined for the multi-processing ringbuffer management. Methods are provided to build the setup dictionaries (necessary parameter objects) for the `Reader`, `Writer` or `Observer` instances, respectively. Further, methods are provided to allow an index processing (e.g. listeners) and to pause data processing.

Index processing overview: to achieve a proper handling for writing and reading data into a ringbuffer slot regarding the possible definition of more than one reader for a ringbuffer (multiprocessing mode on function level) the determination of the correct slot index is essential. For a `Writer` or `Reader` two `SimpleQueues` are defined, respectively. In principle one queue contains the index to be used next and the other the index just processed. An `Observer` is treated in a simpler way by using an own global index variable for each instance.

The queues are passed via the setup dictionary to the corresponding `Reader` or `Writer` instance, respectively, and processed there. In addition, the buffer manager provides methods to control the index determination (`_writer/_reader/_observer_queue_listener()`), started in own threads.

Writer index:

- `writer_empty_queue`: contains the slot numbers (initially filled); defines the next free ringbuffer slot.
  - fetched (removed) in the class `Writer` -> `get_new_buffer()`

- last processed slot number is refilled in `_increment_reader_pointer()`
- `writer_filled_queue`: empty; contains the last processed slot number (distributed to all defined readers).
  - `process_buffer()`

Reader index:

- manually incremented in `_increment_reader_pointer()` via global variable `read_pointer`
- `done_queue`: empty; already processed slot number
- fetched in `_reader_queue_listener()`
- filled in the class `Reader` -> `get()` via the global variable `_last_get_index`
- `todo_queue`: empty; slot number to be processed next
- fetched in the class `Reader` -> `get()`
- filled in `_writer_queue_listener` (within the list `reader_todo_queue_list`)

Observer index:

- the global variable `obs_pointer` is used; it is an early copy of the `write_pointer` variable
  - defined in `_writer_queue_listener()`
  - directly used as index in `_observeQ_listener()`

important methods:

- `__init__()` constructor to create a new 'FIFO' ringbuffer
- `new_writer()` create new writer
- `new_reader_group()` create reader group
- `new_observer()` create observer
- `buffer_status()` display status: event count, processing rate, occupied slots
- `pause()` disable writer(s) to ringbuffer
- `resume()` (re-)enable writers
- `set_ending()` stop data-taking (gives processes time to finish before shutdown)
- `close()` release shared memory
- `shutdown()` end connected processes, delete ringbuffer

#### **buffer\_status()**

Processing Rate and approximate number of free slots in this ringbuffer. This method is meant for user information purposes only, as the result may not be completely accurate due to race conditions.

##### **Returns**

cumulative event count, number of free slots, processing rate, average deadtime

##### **Return type**

tuple

**new\_observer()**

Method to create a new (Queue based) observer.

Method: a copy of the most recent data (latest write\_pointer) is transferred via a Queue whenever the Queue (of size 1) is empty. Sending data through the Queue is handled in a sparate thread

**Returns**

The `setup_dict` object passed to an `Observer`-instance to give access to the data Queue defined for this ringbuffer.

**Return type**

dict

**new\_reader\_group()**

Method to create a new reader group. The processing workload of a group can be distributed to multiple processes by using the same setup dictionary (`setup_dict`) defined for a `Reader`-object. Each ringbuffer element is processed by one reader group process. It's possible to create multiple reader groups per ringbuffer, where each reader group gets every element written to the ringbuffer. If a reader group is created, at least one `Reader`-class instance MUST steadily call its `get()` method to prevent the ringbuffer from blocking and to allow a safe shutdown.

**Returns**

The `setup_dict` object passed to a `Reader`-instance to grant read access to this ringbuffer.

**Return type**

dict

**new\_writer()**

Method to create a new writer. It is possible to create multiple writers and simply share a setup dictionary definition between different `Writer`-instances (analogues to the behavior of the `new_reader_group`).

**Returns**

The `setup_dict` object passed to a `Writer`-instance to grant write access to this ringbuffer.

**Return type**

dict

**pause()**

Disable writing to ringbuffer (paused)

**resume()**

(Re)enable writing to ringbuffer (resume)

**set\_ending()**

Stop data flow (before shut-down)

**shutdown()**

Shut down the ringbuffer(s): close background threads, terminate associated processes and release the shared memory definitions.

Affect processes using a `Reader`, `Writer` or `Observer` instance to a ringbuffer.

A 'trickle down' approach is used to have as few ringbuffer elements as possible unprocessed. This may not work correctly with more complex signal analysis chains. So always make sure to shut down the ringbuffers in data flow order (start with first element of the chain, the ringbuffer closest to the signal source).

**CAUTION!** If there are loops in the signal analysis chain, this method may end in an infinite loop!

**class** mimocorb.mimo\_buffer.**Observer**(*setup\_dict*)

Class for reading selected elements from a ringbuffer via q multiprocessing Queue

The data transfer is implemented via a multiprocessing Queue and interfaces with the ringbuffer manager (NewBuffer-class).

**get()**

Get latest element from buffer: metadata and data

As new data is provided and transferred as soon as data is read from the Queue, the get() method must not be called too frequently

**class** mimocorb.mimo\_buffer.**Reader**(*setup\_dict*)

Class to read elements from a ringbuffer (multiple-out part).

Ringbuffer elements are structured NumPy arrays and strictly **read-only**. The returned array won't change until the next Reader.get() call is performed, blocking the ringbuffer element for the time being. A program design processing the ringbuffer content has to call the Reader.get()-method in a way that minimizes the ringbuffer lock time.

methods:

- get()
- get\_metadata():

**data\_available()**

Method to check for new data and avoid blocking of consumers

**get()**

Get a new element from the ringbuffer. The last element obtained by calling this function is marked as "processing is done". No memory views of old elements may be accessed after calling this function (memory might change, be corrupted or be inconsistent). This function blocks if there are no new elements in the ringbuffer.

**Raises**

**SystemExit** – When the shutdown()-method of the NewBuffer object has been called, a SystemExit is raised which terminates the process.

**Returns**

One element (structured numpy.ndarray) of the ringbuffer as specified in the NewBuffer()-dtype-object.

**Return type**

numpy.ndarray

**get\_metadata()**

Get the metadata defined for a ringbuffer element of the Reader.get()-method.

**Returns**

Currently a 3-tuple is returned with (counter, timestamp, deadline) which is assigned to the latest element of the ringbuffer. The content of these variables is filled by the Writer-process. The current convention is:

- counter (int): a unique, 0 based, consecutive integer referencing this element
- timestamp (float): the UTC timestamp
- **deadline (float): In a live-data environment, the dead time of the first**  
writer in the analyses chain. This is meant to be the fraction of dead time to active data

capturing time (so 0.0 = no dead time whatsoever; 0.99 = only 1% of the time between this and the last element was spent with active data capturing)

**Return type**

tuple

**class** mimocorb.mimo\_buffer.**Writer**(*setup\_dict*)

Class to write elements into a ringbuffer (multiple-in part).

Ringbuffer elements are structured NumPy arrays. Writing is triggered by a call of `Writer.process_buffer()` or at the next call of `Writer.get_new_buffer()`. The ringbuffer element is blocked while writes to the NumPy array are permitted. A program design processing the ringbuffer content has to call the `Writer.process_buffer()` or `Writer.get_new_buffer()`-methods in a way that minimizes the ringbuffer lock time.

methods:

- `get_new_buffer()`
- `set_metadata()`
- `process_buffer()`

**get\_new\_buffer()**

**Get a new free element in the ringbuffer.**

The last element obtained by calling this function is marked as “ready to be processed”. No memory views of old elements may be accessed after calling this function. This function blocks if there are no free elements in the ringbuffer and always returns a valid NumPy array that can be written to.

**Raises**

**SystemExit** – When the `shutdown()`-method of the `NewBuffer` object has been called, a `SystemExit` is raised which terminates the process.

**Returns**

One free ringbuffer element (structured `numpy.ndarray`) as specified in the `NewBuffer()`-dtype-object. Free elements may contain older data, but they can be safely overwritten.

**Return type**

`numpy.ndarray`

**process\_buffer()**

Mark the current ringbuffer element as “ready to be processed”.

The content of the array **MUST NOT** be changed after calling this function. If there is no current element, nothing happens. As the ringbuffer element is blocked while writing to the NumPy array it is recommended to call `Writer.process_buffer()` as soon as possible to minimize the ringbuffer lock time.

**set\_metadata**(*counter*, *timestamp*, *deadtime*)

Set the metadata defined for the current ringbuffer element. If there is no current ringbuffer element (e.g. because `process_buffer()` has been called or `get_new_buffer()` has not been called yet), nothing happens. Copying metadata from a `Reader` to a `Writer` object (here called `source` and `sink`) can be done with:

```
sink.set_metadata(*source.get_metadata())
```

**Parameters**

- **counter** (*integer* (`np.longlong`)) – a unique, 0 based, consecutive integer referencing this element

- **timestamp** (*float (np.float64)*) – the UTC timestamp
- **deadtime** (*float (np.float64)*) – In a live-data environment, the dead time of the first writer in the analyses chain. This is meant to be the fraction of dead time to active data capturing time (so 0.0 = no dead time whatsoever; 0.99 = only 1% of the time between this and the last element was spent with active data capturing)

Collection of classes to set-up, manage and access ringbuffers and associated functions

```
class mimocorb.buffer_control.buffer_control(buffers_dict, functions_dict, function_config_dict,  
                                             output_directory)
```

Set-up and management ringbuffers and associated sub-processes

Class methods:

- **setup\_buffers()**
- **setup\_workers()**
- **start\_workers()**
- **pause()**
- **resume()**
- **shutdown()**

**display\_functions()**

Print list of functions and buffer associations

**display\_layout()**

Print list of buffers

**pause()**

Pause data acquisition

**resume()**

Re-enable data acquisition after pause

**setup\_workers()**

Set up all the (parallel) worker functions

**shutdown()**

Delete buffers, stop processes by calling the shutdown()-Method of the buffer manager

**start\_workers()**

start all of the (parallel) worker functions

**stop()**

stop writing and reading data, allow processes to finish

```
class mimocorb.buffer_control.rbDrain(source_list=None, config_dict=None, **rb_info)
```

read data from ring buffer and sent to null

```
class mimocorb.buffer_control.rbExport(source_list=None, config_dict=None, **rb_info)
```

Read data from buffer and send to requesting client (via Python yield()). Data are provided by a generator function yielding data and metadata in the `__call__()` method of the class.

```
class mimocorb.buffer_control.rbImport(sink_list=None, config_dict=None, ufunc=None, **rb_info)
```

Read data from external source (e.g. front-end device, file, simulation, etc.) and put data in `mimo_buffer`. Data is read by calling a user-supplied generator function for data and metadata.

**class** mimocorb.buffer\_control.**rbObserver**(*observe\_list=None, config\_dict=None, \*\*rb\_info*)

Deliver data from buffer to an observer process. A tuple (data, metadata) is provided by a generator function ( i.e. via yield()) implemented in the `__call__()` method of the class.

**class** mimocorb.buffer\_control.**rbPut**(*sink\_list=None, config\_dict=None, ufunc=None, \*\*rb\_info*)

Recieve data from external source (e.g. front-end device, file, simulation, etc.) and put data in mimo\_buffer.

Returns False if sink is not active

**class** mimocorb.buffer\_control.**rbTransfer**(*source\_list=None, sink\_list=None, config\_dict=None, ufunc=None, \*\*rb\_info*)

Read data from input buffer, filter data and write to output buffer(s) Data is provided as the argument to a user-defined filter function returning None if data is to be rejected, a number if data is to be copied to another buffer, or a list of processed input data write to additional buffers.

Args:

- buffer configurations (only one source and several sinks, no observers!)
- function `ufunc()` must return
  - None if data to be rejected,
  - int if only raw data to be copied to sink[0]
  - list of parameterized data to be copied to sinks[]

Action:

store accepted data in buffers

**class** mimocorb.buffer\_control.**rbWSObserver**(*observe\_list=None, config\_dict=None, \*\*rb\_info*)

Deliver data from buffer to an observer process

**class** mimocorb.buffer\_control.**rb\_toParquetfile**(*source\_list=None, config\_dict=None, \*\*rb\_info*)

Save data a set of parquet-files packed as a tar archive

**class** mimocorb.buffer\_control.**rb\_toTxtfile**(*source\_list=None, config\_dict=None, \*\*rb\_info*)

Save data to file in csv-format

**class** mimocorb.buffer\_control.**run\_mimoDAQ**(*setup\_filename, verbose=2, debug=False*)

Setup and run Data Aquisition suite with mimoCoRB buffer manager

The layout of ringbuffers and associated functions is defined in a configuration file in yaml format.

Functions:

- `setup()`
- `run()`
- `end()`

Data acquisition stops when either of the following conditions is met:

- number of requested events processed
- requested run-time reached
- inuput source exhausted
- end command issued from Keyboard or graphical interface

**end**(*twait=3.0*)

clean shutdown of daq suite

Arg:

*twait*: waiting time for processes to finish before shutdown

**keyboard\_input**(*cmd\_queue*)

Read keyboard input and send to Queue, running as background-thread to avoid blocking

**class tc**

define terminal color codes

*bufferInfoGUI* Graphical display of buffer status

code adapted from <https://github.com/GuenterQuast/picoDAQ>

`mimocorb.bufferinfoGUI.bufferinfoGUI(Qcmd, Qlog, Qinfo, RBnames=['RB_1'], maxRate=100.0, interval=1000.0)`

Show Buffer Manager logging messages and rate history and command buttons

#### Parameters

- **Qcmd** – multiprocessing.Queue for command passing to calling process
- **Qlog** – multiprocessing.Queue() for logging-info
- **Qinfo** – multiprocessing.Queue() for status info
- **RBnames** – list, buffer names, used as line labels
- **maxrate** – maximum rate for y-axis
- **interval** – update interval for graphics in ms

**class** `mimocorb.bufferinfoGUI.plot_bufferinfo(Q, RBnames, maxRate=1500.0, interval=1000.0)`

display statistics from Buffer Manager

uses multiprocessing.Queue() to display buffer information: total number of events, data acquisition rate, buffer filling level

#### Parameters

- **Q** – multiprocessing.Queue() for status info
- **RBnames** – list, buffer names, used as line labels
- **maxrate** – maximum rate for y-axis
- **interval** – graphics update interval in ms

**plot\_buffer** Collection of classes with graphics functions to plot buffer data

The class `animWaveFormPlotter` is used by the class `plotBuffer()`. The `__call__()` method of this latter class is the entry point to the package.

**class** `mimocorb.plot_buffer.animWaveformPlotter(conf_dict=None, dtypes=None, source_dict=None, fig=None)`

Oscilloscope-like display of wave from buffer data

The `__call__` method of this class receives input data and updates and redraws only the new elements of the figure created in `__init__`

**init()**

plot initial line objects to be animated



```
class mimocorb.plot_buffer.plot_buffer(source_list=None, sink_list=None, observe_list=None,
                                       config_dict=None, **rb_info)
```

Plot data using a mimiCoRB Observer

**histogram\_buffer** collection of classes to produce histograms

Show animated histogram(s) of scalar buffer variable(s)

Because this process runs as a 'Reader' process, the plotting function is executed as a background task in order to avoid blocking of the main task.

The entry point is either the `__call__()` function of class `histogram_buffer`, which connects to a `mimocorb` buffer via the `rbExport` class, or directly the function `plot_Histograms()`, which receives data to be histogrammed via a multiprocessing `Queue()`.

code adapted from <https://github.com/GuenterQuast/picoDAQ>

```
class mimocorb.histogram_buffer.animHists(Hdescr, name='Histograms', fig=None)
```

display histograms, as normalised frequency distributions

```
class mimocorb.histogram_buffer.histogram_buffer(source_list=None, sink_list=None,
                                                  observe_list=None, config_dict=None, **rb_info)
```

Produce Histograms of (scalar) variables in buffer.

Read data from mimiCoRB buffer using the interface class `mimo_control.rbExport` and show histograms of scalar variables selected in the configuration dictionary

Plotting is done by means of the class `plot_Histograms()` running as background process

#### Parameters

- **input** – configuration dictionaries
- **config\_dict** – must contain a block with name 'histograms', formatted as  
                   <name>: [<min>, <max>, <nbins>, <ymin>, <ymax>, <label>, <0/1 for lin/log>]

```
mimocorb.histogram_buffer.plot_Histograms(Q, Hdescriptors, interval, name='Histograms')
```

show animated histogram(s)

Args:

- Q: multiprocessing.Queue()
- Hdescriptors: list of histogram descriptors, where each descriptor is a list:
  - min: minimum value
  - max: maximum value
  - nbins: nubmer of bins
  - ymax: scale factor for highest bin ( $1. = 1/N_{bins}$ )
  - name: name of the quantity being histogrammed
  - type: 0 linear, 1 for logarithmic y scale
- interval: time (in s) between updates
- name: name of histogram window

```
class mimocorb.pulseSimulator.pulseSimulator(config_dict)
```

generate wavevorn data of typical pulses from particle detectors, characterised by an exponential shape with parameters height and mean length tau and a noise contribution.

**init**(*number\_of\_channels=None, number\_of\_values=None, channel\_names=None*)

set parameters from buffer configuration and initialize

**class** mimocorb.parquetReader.**parquetReader**(*config\_dict=None*)

read wavefoem data from parquet file

**init**(*number\_of\_channels=None, number\_of\_values=None, channel\_names=None*)

set parameters from buffer configuration and initialize

**rb\_unittest**: application example for mimo\_buffer

This code may serve as a very basic starting point for own projects

Set-up: 2 ring buffers are defined:

- input Buffer RB\_1: 10 ch x 1024 slots (int32)
- output Buffer RB\_2: 10 ch x 2 slots/ch (float64)

**external\_source**: Template for data import in a mimoCoRB buffer from an external source

Input data is provided as a numpy-array of shape (number\_of\_channels, number\_of\_samples).

**simul\_source.simulation\_source**(*source\_list=None, sink\_list=None, observe\_list=None, config\_dict=None, \*\*rb\_info*)

General example for data import from external source (here: generation of simulated data with module pulseSimulator)

Uses class mimocorb.buffer\_control/rbImport to interface to the newBuffer and Writer classes of the package mimoCoRB.mimo\_buffer

mimiCoRB interacts with this code via a generator (*yield\_data()*), which itself received data via the *\_\_call\_\_* function of the class *dataSource* providing the input data. Configuration parametes in the dictionary *config\_dict* are passed to this class during initialistation. Parameters of the configured buffers are set after after initialisation.

This example may serve as a template for other data sources

### Module **lifetime\_filter**

This (rather complex) module filters waveform data to search for valid signal pulses in the channel data. The goal is to clearly identify coincidences of signals in different layers (indiating the passage of a cosmic ray particle, a muon) and find double-pulse signatures that a muon was stopped in or near a detection layer where the resulting decay-electron produced a delayed pulse. The time difference between the initial and the delayed pulses is the individual lifetime of the muon.

Wave forms passing this filter-criterion an passed on to a new buffer; the decay time and the properties of the signal pulses (height, integral and postition in time) are written to another buffer.

The relevant configuration parameters can be found in the section *calculate\_decay\_time*: of the configuration file.

**lifetime\_filter.calculate\_decay\_time**(*source\_list=None, sink\_list=None, observe\_list=None, config\_dict=None, \*\*rb\_info*)

Calculate decay time as time between double pulses

#### **Input:**

pulse wave forms

#### **Returns:**

None if failed, int or list of pulse parameters if successful

Note: output produced when filter is passed depends on number of defined sinks:

- one sink: input data

- two sinks: input data and double-pulse parameters
- three sinks: input data and double-pulse parameters separately for upwards and for downwards going decay electrons

**plot**: plotting waveforms from buffer using mimoCoRB.buffer\_control.OberserverData

`plot_waveform.plot_waveform(source_list=None, sink_list=None, observe_list=None, config_dict=None, **rb_info)`

Plot waveform data from mimiCoRB buffer

**Parameters**

**input** – configuration dictionary

- **plot\_title**: graphics title to be shown on graph
- **min\_sleeptime**: time between updates
- **sample\_time\_ns**, **channel\_range**, **pretrigger\_samples** and **analogue\_offset** describe the waveform data as for oscilloscope setup

**plot\_histograms**: histogram variable(s) from buffer using mimoCoRB.histogram\_buffer

`plot_histograms.plot_histograms(source_list=None, sink_list=None, observe_list=None, config_dict=None, **rb_info)`

Online display of histogram(s) of variable(s) from mimiCoRB buffer

**Parameters**

**input** – configuration dictionary



## PYTHON MODULE INDEX

### I

`lifetime_filter`, 30

### m

`mimocorb`, 21

`mimocorb.buffer_control`, 26

`mimocorb.bufferinfoGUI`, 28

`mimocorb.histogram_buffer`, 29

`mimocorb.mimo_buffer`, 21

`mimocorb.parquetReader`, 30

`mimocorb.plot_buffer`, 28

`mimocorb.pulseSimulator`, 29

### p

`plot_histograms`, 31

`plot_waveform`, 31

### r

`rb_unittest`, 30

### s

`simul_source`, 30



## A

`animHists` (class in `mimocorb.histogram_buffer`), 29  
`animWaveformPlotter` (class in `mimocorb.plot_buffer`), 28

## B

`buffer_control` (class in `mimocorb.buffer_control`), 26  
`buffer_status()` (`mimocorb.mimo_buffer.NewBuffer` method), 22  
`bufferinfoGUI()` (in module `mimocorb.bufferinfoGUI`), 28

## C

`calculate_decay_time()` (in module `lifetime_filter`), 30

## D

`data_available()` (`mimocorb.mimo_buffer.Reader` method), 24  
`display_functions()` (`mimocorb.buffer_control.buffer_control` method), 26  
`display_layout()` (`mimocorb.buffer_control.buffer_control` method), 26

## E

`end()` (`mimocorb.buffer_control.run_mimoDAQ` method), 27

## G

`get()` (`mimocorb.mimo_buffer.Observer` method), 24  
`get()` (`mimocorb.mimo_buffer.Reader` method), 24  
`get_metadata()` (`mimocorb.mimo_buffer.Reader` method), 24  
`get_new_buffer()` (`mimocorb.mimo_buffer.Writer` method), 25

## H

`histogram_buffer` (class in `mimocorb.histogram_buffer`), 29

## I

`init()` (`mimocorb.parquetReader.parquetReader` method), 30  
`init()` (`mimocorb.plot_buffer.animWaveformPlotter` method), 28  
`init()` (`mimocorb.pulseSimulator.pulseSimulator` method), 29

## K

`keyboard_input()` (`mimocorb.buffer_control.run_mimoDAQ` method), 28

## L

`lifetime_filter`  
 module, 30

## M

`mimocorb`  
 module, 21  
`mimocorb.buffer_control`  
 module, 26  
`mimocorb.bufferinfoGUI`  
 module, 28  
`mimocorb.histogram_buffer`  
 module, 29  
`mimocorb.mimo_buffer`  
 module, 21  
`mimocorb.parquetReader`  
 module, 30  
`mimocorb.plot_buffer`  
 module, 28  
`mimocorb.pulseSimulator`  
 module, 29  
 module  
   `lifetime_filter`, 30  
   `mimocorb`, 21  
   `mimocorb.buffer_control`, 26  
   `mimocorb.bufferinfoGUI`, 28  
   `mimocorb.histogram_buffer`, 29  
   `mimocorb.mimo_buffer`, 21  
   `mimocorb.parquetReader`, 30

mimocorb.plot\_buffer, 28  
mimocorb.pulseSimulator, 29  
plot\_histograms, 31  
plot\_waveform, 31  
rb\_unittest, 30  
simul\_source, 30

## N

new\_observer() (mimocorb.mimo\_buffer.NewBuffer  
method), 22  
new\_reader\_group() (mi-  
mocorb.mimo\_buffer.NewBuffer method),  
23  
new\_writer() (mimocorb.mimo\_buffer.NewBuffer  
method), 23  
NewBuffer (class in mimocorb.mimo\_buffer), 21

## O

Observer (class in mimocorb.mimo\_buffer), 23

## P

parquetReader (class in mimocorb.parquetReader), 30  
pause() (mimocorb.buffer\_control.buffer\_control  
method), 26  
pause() (mimocorb.mimo\_buffer.NewBuffer method), 23  
plot\_buffer (class in mimocorb.plot\_buffer), 28  
plot\_bufferinfo (class in mimocorb.bufferinfoGUI),  
28  
plot\_histograms  
module, 31  
plot\_Histograms() (in module mi-  
mocorb.histogram\_buffer), 29  
plot\_histograms() (in module plot\_histograms), 31  
plot\_waveform  
module, 31  
plot\_waveform() (in module plot\_waveform), 31  
process\_buffer() (mimocorb.mimo\_buffer.Writer  
method), 25  
pulseSimulator (class in mimocorb.pulseSimulator),  
29

## R

rb\_toParquetfile (class in mimocorb.buffer\_control),  
27  
rb\_toTxtfile (class in mimocorb.buffer\_control), 27  
rb\_unittest  
module, 30  
rbDrain (class in mimocorb.buffer\_control), 26  
rbExport (class in mimocorb.buffer\_control), 26  
rbImport (class in mimocorb.buffer\_control), 26  
rbObserver (class in mimocorb.buffer\_control), 26  
rbPut (class in mimocorb.buffer\_control), 27  
rbTransfer (class in mimocorb.buffer\_control), 27  
rbWSObserver (class in mimocorb.buffer\_control), 27

Reader (class in mimocorb.mimo\_buffer), 24  
resume() (mimocorb.buffer\_control.buffer\_control  
method), 26  
resume() (mimocorb.mimo\_buffer.NewBuffer method),  
23  
run\_mimoDAQ (class in mimocorb.buffer\_control), 27  
run\_mimoDAQ.tc (class in mimocorb.buffer\_control), 28

## S

set\_ending() (mimocorb.mimo\_buffer.NewBuffer  
method), 23  
set\_metadata() (mimocorb.mimo\_buffer.Writer  
method), 25  
setup\_workers() (mi-  
mocorb.buffer\_control.buffer\_control method),  
26  
shutdown() (mimocorb.buffer\_control.buffer\_control  
method), 26  
shutdown() (mimocorb.mimo\_buffer.NewBuffer  
method), 23  
simul\_source  
module, 30  
simulation\_source() (in module simul\_source), 30  
start\_workers() (mi-  
mocorb.buffer\_control.buffer\_control method),  
26  
stop() (mimocorb.buffer\_control.buffer\_control  
method), 26

## W

Writer (class in mimocorb.mimo\_buffer), 25